

LENOXINS

Manufacturers of Fine Quality Dental Instruments

**NEW
EDITION**

SURGICAL

HYGIENISTS

ORTHODONTIC

DESIGN
PRECISION
QUALITY

For Orders

+ 92 - 52 - 3241123

+ 92 - 320 - 6338182

lenoxsurgical@gmail.com

www.lenoxsurgical.com



Best Selling
Dental Instruments
Catalogue



LENOX SURGICAL

A Little About Us:

LENOX Surgical ins. Present This New Edition Dental Instruments Catalogue With Hope That It Would Be An Excellent Guide To All Over Importer, Distributors, And Dental Professionals Worldwide.

Lenox Surgical Instruments is a ISO-9001-2015, ISO-13485 -2016CE, cGMP (FDA-USA) Certified Company and Head Office is Located in Industrial Area of Sialkot, Which is Well Known Worldwide For Manufacturing Quality Surgical & Dental Instruments.

We Are Also Very Grateful for The Loyal Support Given By Our Customer/ Distributors, Who Have Been Selling Our Instruments in Their Market With Pride And Honor.

We Are Pledged To Serve More In The Years Ahead To Our Quality Minded Customers With Much More Higher Standard of Instruments and Services.

We Offer A Comprehensive Portfolio of All Dental Specialities, with Ready Stock for Quick Deliveries.

We Also Welcome Your Enquires for Instruments Not Listed in This Catalogue Both in Surgical & Dental.

At LENOX We Feel That Good Communication is Key To Any Successful Business Relationship, Please Feel Free To Contact Us By The Methods That Are The Most Convenient For Your Business And We Will Do Our Best To Always Provide Great Quality Products And Services.

Your Sincerely,
LENOX SURGICAL INSTRUMENTS
Sialkot, Pakistan.



LENOX SUR warranties all its product against poor craftsmanship and materials. Any product found defective will be repaired or replaced free of charge. This warranty does not cover any products which become defective through misuse, negligence and mishandling.

LENOX SURGICAL INSTRUMENTS



Index



01	Extraction Forceps (English American Children)		Page 01-09
02	Bone Rongeurs		Page 10
03	Root Elevators		Page 11-16
04	Dental Anesthesia Syringes		Page 17
05	Handles & Mouth Mirrors		Page 18
06	Dental Tweezers		Page 19
07	Dental Scissors		Page 20
08	Micro Surgery / Needle Holders		Page 21-23
09	Explorers & Probes		Page 25-26
10	Excavators / Pluggers / Condensers / Carvers		Page 27
11	Burnishers / Hollenback / Carving Knives		Page 28
12	Filling / Composite		Page 29-30
13	Dental Scalars & Curettes		Page 31-33
14	Periosteal Elevators		Page 34
15	Dental General Instruments		Page 35
16	Rubber Dam Instruments		Page 36
17	Dental Crown Instruments		Page 37
18	Dental Implants Instruments		Page 38-40
19	Orthodontics Instruments		Page 41-43
20	Impression Trays		Page 44
21	Dental Hollow Wares / Cassettes		Page 45-46




ORAL SURGERY

Extraction Forceps

CHILDREN & ADULT TEETH

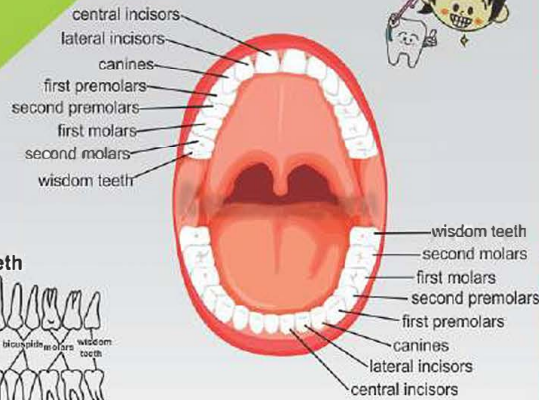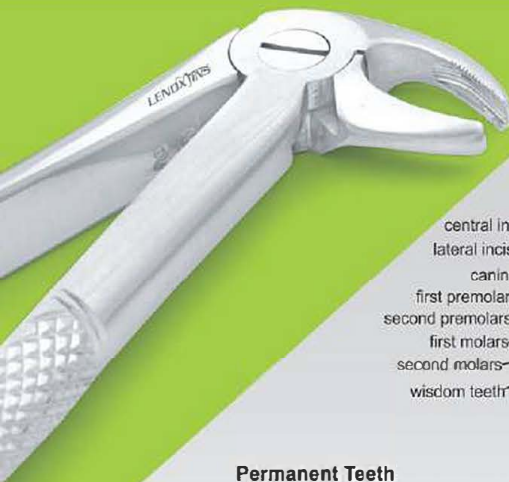

Permanent Teeth

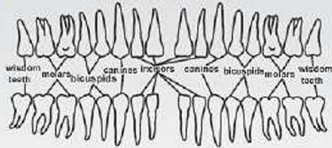

Requirement of the Dental Forceps INTRA-ALVEOLAR EXTRACTION

Dental extraction forceps has three components which is blades and handles united together by a hinge joints.

1-The serrated handle must possess a suitable size to rest comfortably in the operator's hand and should long enough to afford use to strong and steady extraction movements. Also it must be suitable in shape and designs to suit the area of the tooth.



AMERICAN PATTERN HANDLE



ENGLISH PATTERN HANDLE



PERFORATED HANDLE



ANATOMICAL HANDLE

01

Dental Extraction Set



- 01 LD13-104- Cotton Plier
- 02 LD1 I-108- Mirror Handle Cone Socket
- 03 LD17-111- Minnesota Check Retractor
- 04 LD2-108- Extracting Forceps Lower Anterior
- 05 LD1-125- Extracting Forceps Lower Molars
- 06 LD4-108- Extracting Forceps Upper Molars
- 07 LD2-103- Extracting Forceps Upper Pre-Molar
- 08 LD7-103- Anesthesia Syringe Self Aspirating
- 09 LD17-117- Austin Tissue Retractor
- 10 LD15-102- Crile Wood Needle Holder 15cm
- 11 LD16- Pean Hemostat Forceps Long (Cvd 16cm)
- 12 LD6-136- Elevator Coupland 4mm
- 13 LD6-113- Elevator Apexo Str 3mm
- 14 LD6-129- Elevator Cryer L
- 15 LD6-130- Elevator Cryer R
- 16 LD6138- Elevator Crana Pick E8
- 17 LD5-102- Bone Rongeur
- 18 LD12-106 Mirror Front Surface #5 (pack of 12)
- 19 LD19-101- Bone File Cross Cut
- 20 LSP-101- Periosteal Elevator Molt
- 21 LSM-103- Lucas Bone Curette Serrated 3mm
- 22 LSM-121- Lucas Bone Curette Non Serrated 3mm
- 23 LD6-164- Heidbrink Root Tip Pick Double Ended
- 24 LD14-120- Kelly Super Cut Scissors Cvd 16cm
- 25 LD18-106- Scalpel Handle # 3 With Ruler

Dental Extraction Set



- 01 LD13-108- Adson Brown Tissue Forceps 12cm
- 02 LD5-103- Blumenthal Rongeur 15cm
- 03 LDT2-108- Goldman-Iox Tc Scissors Cvd 13cm
- 04 LD15-102- Crile-wood Needle Holder 15cm
- 05 LD7-103- Self Aspirating Syringe 1.8ml
- 06 LD5-301- #301 Apexo Elevator Str 2mm
- 07 LD6-180- Arrow PDL Black
- 08 LD19-101- Bone File
- 09 LPS-120- Molt Periosteal Elevator
- 10 LSM-103- Lucas Curette 3mm Serrated
- 11 LD18-102- Long Scalpel Handle #7 Angled

Extracting Forceps (Set of 6 Pcs)

Cross serrations for a non-slip grasp. Easily enter the sub gingival.

Universal Atraumatic Extraction Forceps Set

Unique Design.



Lower Anterior

LDA1-153GS



Upper Anterior

LDA1-154GS



Upper Premolars

LDA1-155GS



Lower Premolars

LDA1-156GS



Lower Molars

LDA1-157GS



Upper Molars

LDA1-158GS



02

Extracting Forceps (Set of 6 Pcs)

Cross serrations for a non-slip grasp. Easily enter the sub gingival.

Black Anti-Reflection Coated

Universal Atraumatic Extraction Forceps Set

Unique Design.



Lower Anterior

LDA1-153BG



Upper Anterior

LDA1-154BG



Upper Premolars

LDA1-155BG



Lower Premolars

LDA1-156BG



Lower Molars

LDA1-157BG



Upper Molars

LDA1-158BG



03



LD1-102
Fig.2, Upper Laterals
and Canines



LD1-103
Fig.3
Upper Incisors
and Canines



LD1-104
Fig.4
Lower Incisors
and Canines



LD1-106
Fig.8
Lower Premolars



LD1-110
Fig.19
Upper Third
Molars



LD1-111
Fig.20
Lower Third
Molars



LD1-112
Fig.21
Lower Molars



LD1-114
Fig.30
Upper Roots



LD1-115
33-Lower
Roots



LD1-116
33A-Lower
Roots



LD1-117
73A-Lower Molars



LD1-119
45-Lower roots



LD1-120
44-Upper roots



LD1-200
153-Root fragments



LD1-127
74N-Lower Roots



LD1-123
Routurier (# 22 1/2 L)
Lower Molars and
wisdoms, left

**WISDOM
EXTRACTING
FORCEPS**



LD1-124
Routurier (# 22 1/2 R)
Lower Molars and
wisdoms, right



LD1-220
97-Upper Roots
very fine



LD1-125
#73-Lower Molars



LD1-128
#86A-Lower Molars

DEEP GRIPPING FORCEPS

- ▶ Parallel gripping beak
- ▶ Aggressive edge serrations
- ▶ Long axis serrations within the beaks
- ▶ Thick beaks for superior penetration
- ▶ Contact the root remaining coronal structure
- ▶ Used on partially intact, broken down or sectioned teeth



LENOX INSTRUMENTS

Oral Surgery



No Fracture | No Tissue Damage | Extraction with Easy | Slim Light Weight Pattern



**LD1-129
35**
Upper Premolars
Deep-gripping



Thin Tip



**LD1-130
34**
Upper Incisors and
Premolars
Deep-gripping



Wide Tip



**LD1-131
34A**
Upper Centrals / Cains
Deep-gripping



Wide Tip



**LD1-132
35A**
Upper Incisors and
Premolars
Deep-gripping



Thin Tip



**LD1-133
36**
Lower Premolars
Deep-gripping



Wide Tip




**LD1-134
36A**
Lower Premolars
Deep-gripping




**LD1-135
336**
Lower Premolars
Deep-gripping



**LD1-136
F-7C**
Lower Molar

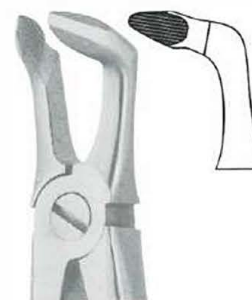

**LD1-137
9**
Apex Forceps
Lower Universal
Deep Root Pick

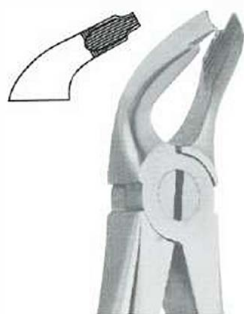


**LD1-138
79N**
Lower Third Molars
Deep-gripping



**LD1-139
79N**
Lower Third Molars



**LD1-140
67**
Upper Third Molars

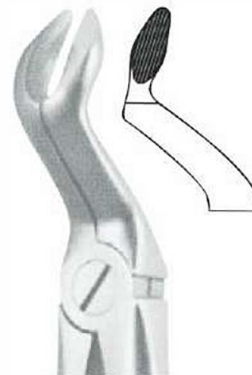

**LD1-141
67A**
Upper Third Molars



**LD1-142
67N**
Upper Third Molars
Deep-gripping

04

Extracting Forceps (Set of 10 Pcs)



05



LD1-143 Fig.1, Upper Laterals and Canines



LD1-144
Fig. 7
Upper Premolars



LD1-145
Fig. 13
lower Premolars



LD1-146
Fig. 17 Upper Molars Right



LD1-147
Fig. 18 Upper Molars Left



LD1-148
Fig. 22
Lower Molars



LD1-149
Fig. 79, Lower Third Molars



LD1-150
Fig. 74
Lower Roots



LD1-151
Fig. 67A
Upper Third Molars



LD1-152
Fig. 51A
Upper Roots



LDS1-10- Set of 10 Pieces
Fig. 1, Fig. 7, Fig. 13, Fig. 17
Fig. 18, Fig. 22, Fig. 79, Fig. 74
Fig. 67A, Fig. 51A

Universal Pediatric Forceps



LD1-153
Fig. 365/150SK
Upper anteriors and Roots



LD1-154
Fig. 360/17SK
Lower Molars



LD1-155
Fig. 370/151SK
Lower anteriors and Roots

Root Pick Extracting Forceps



LD1-156
Forcep Fig. 301



LD1-157
Forcep Fig. 300

Perforated Light Weight Handles
with Strong Serrated Gripping Action.



LD2-101
Fig. 1
Upper Anterior

Side View



LD2-101
Set of 10 Pieces
Fig. 1, Fig. 23, Fig. 74
Fig. 88L, Fig. 88R, Fig. 74
Fig. 106, Fig. 174, Fig. 859

LD2-102
Fig. 23 Cow Horn
Lower Molar

LD2-103
Fig. 35
Upper Premolar

LD2-104
Fig. 88L
Upper Molar Left

LD2-105
Fig. 88R
Upper Molar Right

LD2-106
Fig. 106
Lower Incisors

LD2-107
Fig. 74
Lower Molars

LD2-108
Fig. 174
Extra Narrow
Lower Anterior
Apical Retention
Forceps

LD2-109
Fig. 67A
Upper Molars

LD2-110
Lower Roots

LD2-111
Upper Roots

LD2-112
Lower Roots



BABY EXTRACTING FORCEPS

To avoid fear of small patients, LENOX INC forceps are that small to fit completely inside a hand. For an easy use they are additionally equipped with a spring, LENOX INC extracting forceps assure a firm, accurate grasp of primary tooth and root surface for smooth, seamless extraction reducing the risk of complications. Hand crafted by Surgical Grade Non-Rusty Stainless steel superior performance and longevity.


Think Teeth

Keep up good dental habits at every age.

LENOX *INSTRUMENTS*



Pediatric Tooth Extracting Forceps Sets

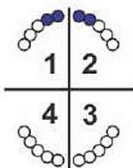

LD3-101
English Pattern
Klein #137 Pedodontic,
Upper Incisors
Spring handle



LD3-102
English Pattern Klein #5
Pedodontic Lower Incisors.



LD3-103
English Pattern Klein #7
Pedodontic Lower Roots,



LD3-104
English Pattern Klein #6
Pedodontic Lower Molars,
Spring Handle



LD3-105
English Pattern Klein #139
Pedodontic Upper Premolars



LD3-106
English Pattern Klein #51S
Pedodontic Upper root




LD3-107
English Pattern Klein #3
Pedodontic Upper
Molars Universal



LDS4-120
Set of 7 Pieces

Pediatric Tooth Extracting Forceps Sets

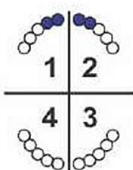

LD3-108
English Pattern
Klein #137 Pedodontic,
Upper Incisors
Spring handle



LD3-109
English Pattern Klein #5
Pedodontic Lower Incisors.



LD3-110
English Pattern Klein #7
Pedodontic Lower Roots,



LD3-111
English Pattern Klein #6
Pedodontic Lower Molars,
Spring Handle



LD3-112
English Pattern Klein #139
Pedodontic Upper Premolars



LD3-113
English Pattern Klein #51S
Pedodontic Upper root



LD3-114
English Pattern Klein #3
Pedodontic Upper
Molars Universal



LDS5-130
Set of 7 Pieces



LD4-101
Fig. 62
incisors, premolars,
deciduous teeth



LD4-102
#10S
upper roots



LD4-103
Harris
Fig. 18R



LD4-104
Harris
Fig. 18L



LD4-105
Fig. 53R
upper molars, right



LD4-106
Fig. 53L
upper molars, left



LD4-107
NEVIUS
Fig. 88R
upper molars, right



LD4-108
Fig. 24
upper molars



LD4-109
#210S
upper third molars



LD4-110
#99C Kells
upper incisors,
canines, premolars



LD4-111
Nevius
Fig. 88L
upper molars, left



LD4-112
Tomes Fig. 69
upper and lower
roots



LD4-113
Hull Fig. 101
Deciduous teeth



LD4-114
Fig. 150 AS



09



LD4-116
Fig. 17
Lower molars



LD4-117
Fig. 16
Harris Lower molars,
COW HORN beak



LD4-118
Fig. 222
Lower wisdoms



LD4-119
Fig. 69
Tomes Fragments,
uppers and lowers



LD4-115
Fig. 222
Lower wisdoms

LD4-120
Fig. 286
Upper roots,
premolars and incisors

LD4-121
Fig. 171
Wisdom and molar roots

LD4-122
Fig. 166
Universal Application

LD4-123
Fig. 169
Centrals laterals premolars
lower canines and all roots



LD4-124
Fig.150 Apical
Upper Universal Forceps



LD4-125
Fig.151 Apical
Lower Universal Forceps



LD4-127
Tweezer Albrecht



LD4-128
Scissor Handle
Lower Wisdoms

BONE RONGEUR

INSTRUMENTS FOR CUTTING BONES

- Most commonly used
- Sharp blades- squeezed together cutting/pinching through bone
- Leaf spring between the handle : instrument opens when hand pressure is released
- Repeated cuts without manually reopening



Oral Surgery

LENOXINS



LD5-101
Blumenthal 90°
15cm



LD5-102
Blumenthal 45°
15cm



LD5-103
Blumenthal 30°
15cm



LD5-104
Micro Friedman
La-grange Bend
14cm



LD5-105
Mini Friedman
14.5cm



LD5-106
Mini Friedman 90°
14.5cm

SOFT TISSUE NIPPERS



LD5-115 - 10cm
Compact, sharp nippers
with smooth action useful
for accentuating interproximal
gingival contour during
soft tissue surgery.

DOUBLE ACTION RONGEURS



LD5-107
Double Action
Side-Cutting
17 cm




Boehler
Double Action
STR
LD5-108 - 15 cm-2mm
LD5-109 - 15 cm-3mm
CVD
LD5-110 - 15 cm-2mm
LD5-111 - 15 cm-3mm



LD5-112
Beyer
Double Action
18 cm

KERISON RONGEUR

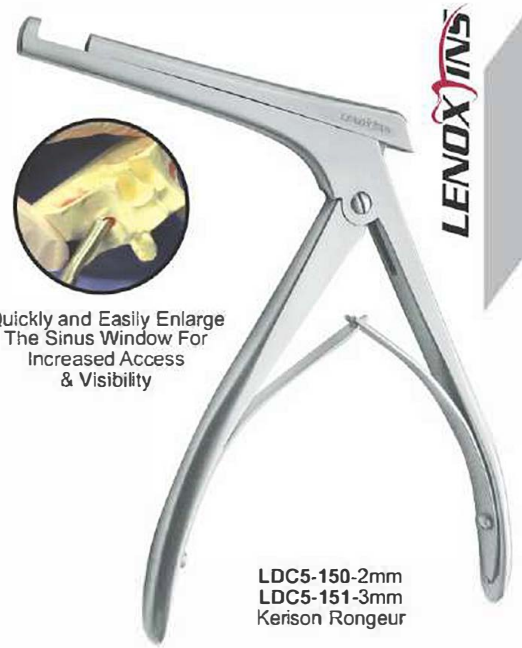


Quickly and Easily Enlarge
The Sinus Window For
Increased Access
& Visibility

LDC5-150-2mm
LDC5-151-3mm
Kerison Rongeur

HANDLES FOR ROOT ELEVATORS


LENOX INS

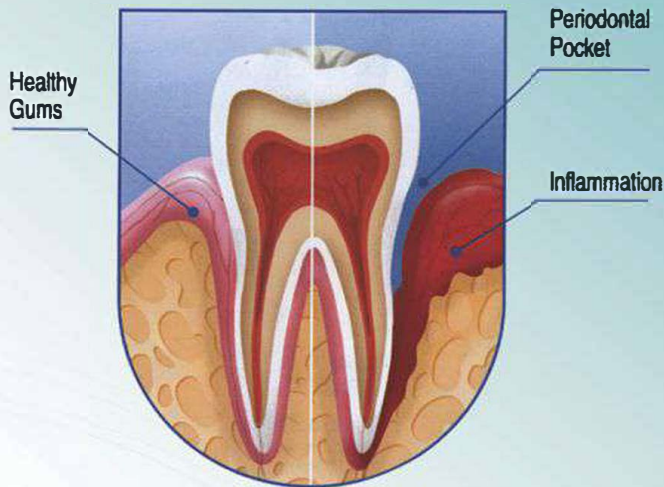

Elevators Classification

1- According to use:

- 1-Elevators designed to remove the entire tooth, straight elevators, hospital pattern and winter elevator
- 2-Elevators designed to remove roots broken off at the gingival line e.g. Apexo elevator, Coupland and lido lavien elevators
- 3-Elevators designed to remove roots broken off half way to the apex e.g. curved elevator hospital pattern, winter elevator and Apexo elevator

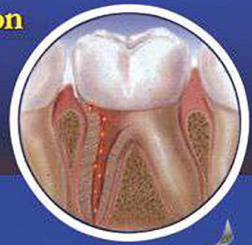

HANDLES FOR ROOT ELEVATORS POINTS / A-B-C-D-E-F-G-H



11

Flexible Periotomes

Gold Titanium Coated Points With Micro Serration



LD6-102
Flexible Periotome 2mm STR



Root Screw Extractor



LENOX INS Flexitome (Set of 6 Pcs)

Super Flexibility, Micro Serrated, Special Alloy, Hard To Be Broken



LD6-105
Straight 2mm



LD6-106
Curved 2mm



LD6-107
Curved 3mm



LD6-108
Side Cutting
2mm



LD6-109
Back Cutting
2mm

Super Flexibility, None Serrated, Special Alloy, Hard To Be Broken

LD6-105NS
Straight 2mm

LD6-106NS
Curved 2mm

LD6-107NS
Curved 3mm

LD6-108NS
Side Cutting
2mm

LD6-109NS
Back Cutting
2mm

LD6-104
Straight 3mm
LD6-104NS
Straight 3mm

LD6-06
Flex-PDL Set of 6




Flexitome Cassette


Apexo Root Elevators



Bein Root Elevators

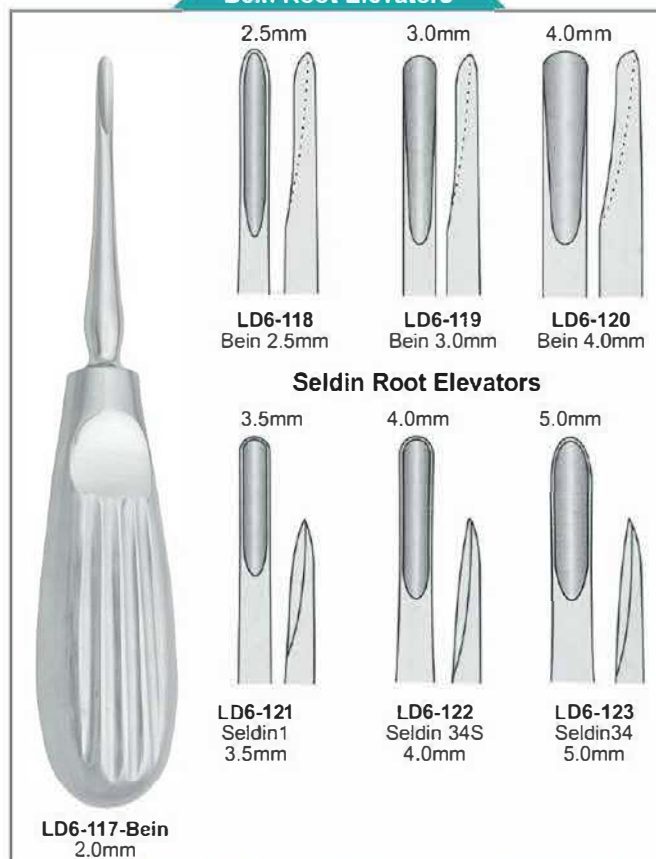

Seldin Root Elevators

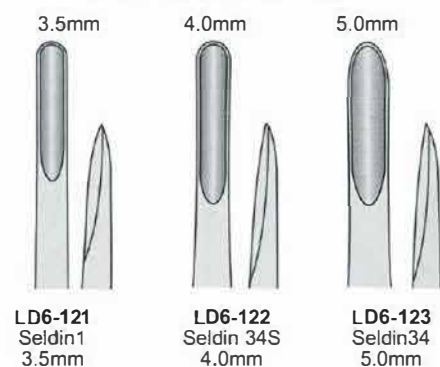

Apical Root Elevators



Cryer Root Elevators



Coupland Root Elevators



Crane Pick



Lindo Root Elevators

Bernard Root Elevators




Claw-Lift Roots Elevator

Heidbrink Root Tip Picks



For those clinicians who are interested the entire luxating Elevators line. **LENOX INS** has created sets with seven micro serrated instruments for atraumatic extraction.

PDL Elevators (Set of 7 Pcs)

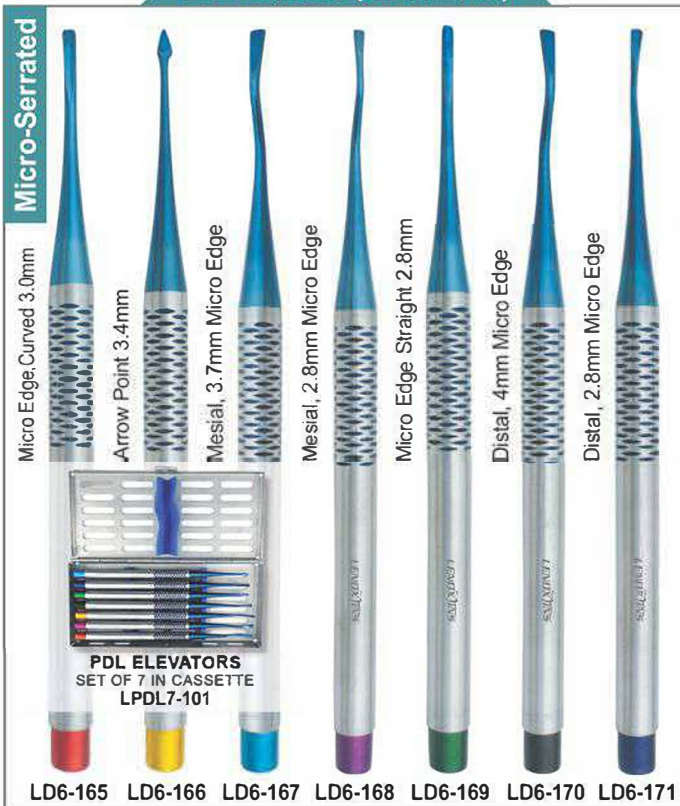

PDL Elevators (Set of 7 Pcs)




PDL Elevators (Set of 7 Pcs)



Bernard Flexible (Set of 4 Pcs)



Luxation Elevators

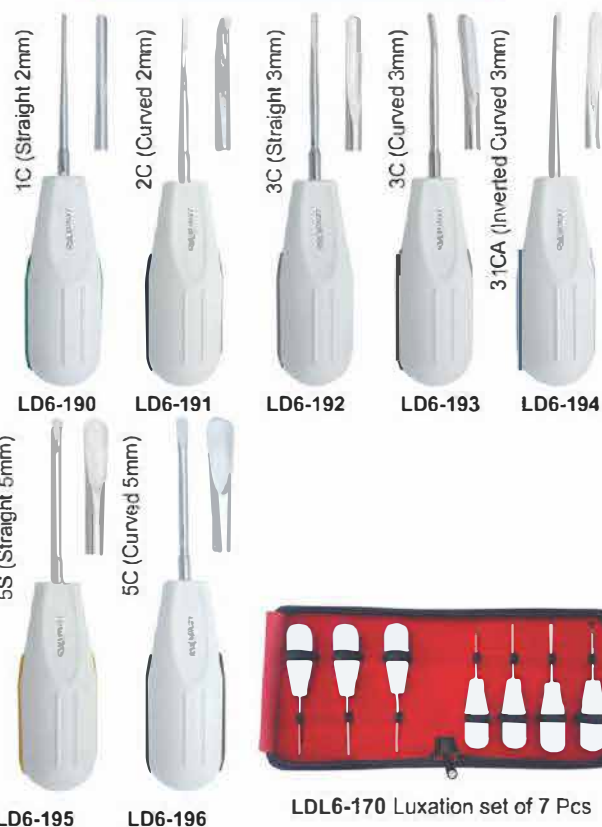

Ergonomic Handles Elevators Set



Potts Luxation Elevators



Chompret Syndesmote



Surgical Curettes



Lucas Surgical Curettes



Laxation Set, Delicated Points



LD6-208
Straight blade
2mm



LD6-209
Curved blade
3mm



LD6-210
Straight blade
3mm



LD6-211
Contra angle
3mm



LD6-212
Inverted curved
3mm



LD6-213
Straight blade
5mm



LD6-214
Curved blade
5mm



LD6-215
Straight dual edge
3mm / 5mm



LD6-216
Straight dual edge
1.5mm / 3mm



LD6-217
Right to bend
3mm



LD6-218
Left to bend
3mm



LD6-219
Right to the
bending 5mm



LD6-220
Left to bend
5mm

Lucas-X curettes



Fig 1, 112.5°, Angle

LXC-1-2mm
LXC-2-3mm
LXC-3-5mm



Fig 2, 135°, Angle

LXC-4-2mm
LXC-5-3mm
LXC-6-5mm



Fig 3, 157.5°, Angle

LXC-7-2mm
LXC-8-3mm
LXC-9-5mm

Micro Serrated

Bone Files



LBF-101
Miller
Cross Cut



LBF-102
Miller 2X
Straight Cut



LBF-103
Wahl 1



LD6-245



LD6-246



LD6-247



LD6-248




LDL6-250

Root Fragment Removal elevators, Set of 4 (Twisted Style)



Self Aspirating Local Anesthesia Syringes

Suitable For Use with All 1.8ml
Dental Cartridges

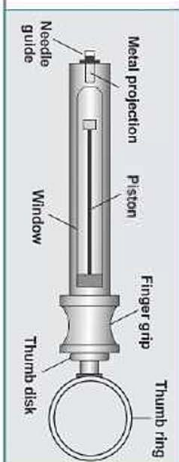

LD7-101



LD7-102




LD7-103



LD7-104




LD7-105

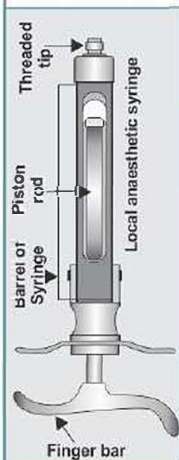

Mechanical Aspirating Local Anesthesia Syringes

Suitable For Use with All 1.8ml
Dental Cartridges

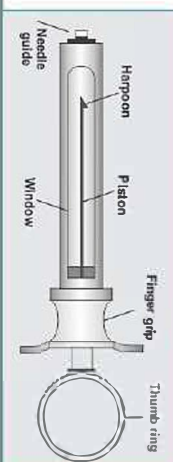

LD8-101



LD8-102




LD8-103



LD8-104



LD8-105

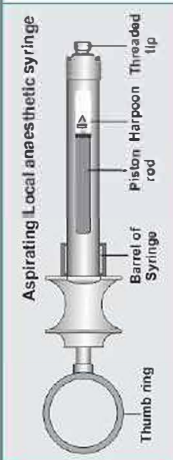

Intraligamental Syringes



LD9-1001



LD9-1002



LD9-1003



Made Of Stainless Steel Color Coated. Suitable For Use with All 1.8ml
Dental Cartridges. Inject 0.05ml with Every Click. Exact Dosage Every Time.
Supplied with Angled Tubes.

Gun Style Intraligamental Syringes

Use Lenox's Premium Quality Stainless Steel Syringe to give
the anesthesia to patients without any problem.



LDG-10 STRAIGHT
LDG-11 ANGLED

Made of Stainless Steel
Lenox Ject
Metric Thread
Inject 0.06ml with every click
Exact dosage every time
Suitable for use with all 1.8ml dental Cartridges
Supplied with Straight and angled tubes.

Handles for Mouth mirrors

LD11-101 Mirror Handle Stem Round (6mm) Length 13cm



LD11-102 Mirror Handle Stem octagonal (6mm) Length 13cm



LD11-103 Mirror Handle Slim Shape Pattern (6mm) Length 15cm



LD11-104 Hollow Mirror Handle Stem Extra Light (7.5mm) Length 13cm



LD11-105 Mirror Handle Stem with Rule 35mm (6mm) Length 13cm



LD11-106 Mirror Handle Stem with Rule 35mm (6mm) Length 13cm (Cone Socket)



LD11-107 Hollow Mirror Handle Stem Extra Light (7.5mm) Length 13cm (Cone Socket)



LD11-108 Hollow Mirror Handle Slim Shape Pattern (9mm) Length 13cm (Cone Socket)



LD11-109 Hollow Mirror Handle Stem Extra Light (9mm) Length 13cm (Cone Socket)



LD11-110 Hollow Mirror Handle blue Titanium Coated (9mm) Length 13cm (Cone Socket)



LD11-111 Hollow Mirror Handle Light weight(9.5mm) Length 13cm (Cone Socket)



Dental Mouth Mirrors Black Colour Coated



LD11-112B-#4 (22mm)
LD11-114B-#5 (24mm)
Combo of Double Ended Black HD Mirror, With Blue Titanium Coated Hollow Handle.

Dental Mouth Mirrors

Double Sided Mouth Mirror

Mouth Mirror plain



LD12-101 Fig. 3 Ø20mm
LD12-102 Fig. 4 Ø22mm
LD12-103 Fig. 5 Ø24mm

Mouth Mirror Cone Socket Front Surface



LD12-104 Fig. 3 Ø20mm
LD12-105 Fig. 4 Ø22mm
LD12-106 Fig. 5 Ø24mm

Front Surface



LD12-107 #4 (22mm)
LD12-108 #5 (24mm)
LD12-109 #4 (22mm)
LD12-110 #5 (24mm)

Dental Tweezers

- Tweezers are used particularly to handle and place cotton pellets, wedges and articulation
- Tweezers are important Dental Lab Equipment, made of non-magnetic stainless steel.
- They are designed particularly to hold the tiniest article, provide better grip on the jaw.



Oral Surgery



LD13-101
London College
15cm



LD13-102
Meriam
16cm



LD13-103
London College
15cm



LD13-104
Meriam
16cm



LD13-105
14cm 1x2
Tissue Forcep



LD13-106
College with Lock
15cm Cvd.



LD13-107
Hemingway
15cm



LD13-108
Adson Brown
12cm



LD13-109
Adson
12cm



LD13-110
Adson
12cm 1x2



LD13-111
Bonding Bracket
Holder, 13.5cm



LD13-112
Bonding Bracket
Holder, 13.5cm



LD13-113
Bonding Bracket
Holder, 13.5cm



LD13-114
Miller
15cm Str.



LD13-115
Miller
15cm Cvd.



LD13-116
Bonding Bracket
Holder, 13.5cm



LD13-117
Bonding Bracket
Holder, 13.5cm



LD13-118
Bonding Bracket
Holder, 13.5cm



LD13-119
Bonding Bracket
Holder, 13.5cm



LD13-120
Com Suture
Pliers, 15cm



LD13-121
Tissue grasping
forceps, atraumatic
17.5 cm



LD13-122-12cm
LD13-123-15cm
Micro Adson



LD13-124
12cm Mini Adson
CVD




LD13-135
Cotton Plier
Wide Grip



LD13-140
16 cm
Cotton Plier
Fine Point

Lab Tweezers

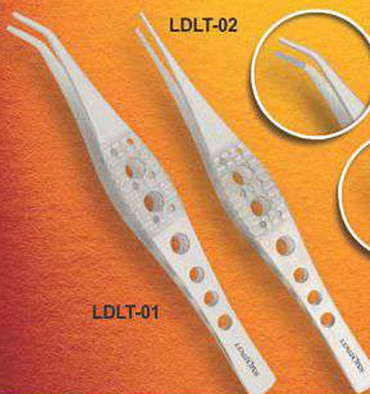

- LDLT-01** Forceps 14cm Curved, Ergonomic Paddle Handle, 1.5mm tip, smooth.
- LDLT-02** Forceps 14cm Straight, Ergonomic Paddle Handle, 1.5mm tip, smooth.
- LDLT-03** Forceps 14cm Curved, Ergonomic Paddle Handle, 1.5mm tip, serrated.
- LDLT-04** Forceps 14cm Straight, Ergonomic Paddle Handle, 1.5mm tip, serrated.
- LDLT-05** Forceps 14cm Straight, Ergonomic Paddle Handle, Needle tip.

DENTAL SCISSORS

Oral Surgery

LENOX INS

Functions and Features:

- Scissors designed for suture removal.
- Blades – one has a curved edge and one has a notched edge.
- During suture removal, suture is scooped up by blade with the notched end before cutting to prevent tissue trauma.

General Scissors



Scissors Tungsten Carbide

Tungsten Carbide scissors feature cutting edges of permanently welded Tungsten carbide inside. These Scissors will retain a fine sharp cutting edge for long time and will outlast any standard stainless steel dental scissors. Both handles are gold-plated for easy identification.



Super Cut Scissors

Micro Serrated



T.C Super Cut Scissors

Tungston carbide inside, One Blade Micro Serrated



Castroviejo Needle Holder
Diamond Dust Coated


LDTTC15-12
CVD 14cm
LDTTC15-13
CVD 16cm
LDTTC15-14
CVD 18cm



LDTTC15-18
CVD 14cm
LDTTC15-19
CVD 16cm
LDTTC15-20
CVD 18cm



Castroviejo TC Twisted Joint
LDTTC15-500 STR 18cm
LDTTC15-510 CVD 18cm
LDTTC15-502 STR 16cm
LDTTC15-512 CVD 16cm



Castroviejo TC Fine Pointed
LDTTC15-400 CVD 18cm
LDTTC15-410 STR 18cm
LDTTC15-402 CVD 16cm
LDTTC15-412 STR 16cm



Micro Needle Holder TC
LDTTC15-425 STR 16cm
LDTTC15-427 CVD 16cm



Castroviejo Fine Pointed
Twisted Joint
LDDDD15-520 STR 18cm
LDDDD15-522 CVD 18cm
LDDDD15-524 STR 16cm
LDDDD15-526 CVD 16cm

Castroviejo TC
LDTTC15-09 STR 14cm
LDTTC15-10 STR 16cm
LDTTC15-11 STR 18cm

Castroviejo TC Fine Pointed
LDTTC15-15 STR 14cm
LDTTC15-16 STR 16cm
LDTTC15-17 STR 18cm

Flex Spring Scissors
Castroviejo Micro Scissors TC


LD14-110 - 15cm
LD14-111 - 17cm
Castroviejo
Scissor, Straight

LD14-112
15cm
LD14-113
17cm



LD14-114
Noyes
11cm
Straight

Curved



LD14-116
Noyes
11cm
Angular



Castroviejo TC Scissor
LDTTC2-101 - 15cm
LDTTC2-102 - 17cm

Curved
LDTTC2-103
15cm
LDTTC2-104
17cm



Castroviejo TC Scissor
LDTTC2-400 - 18cm, STR
LDTTC2-410 - 18cm, CVD
LDTTC2-402 - 16cm, STR
LDTTC2-412 - 16cm, CVD



LDTTC2-450
Micro Scissor TC
90° Curved, 18 cm



LDTTC2-455
Micro Scissor TC
90° Angled, 18 cm

Suture Scissor Duck Blade



Suture Scissor Duck Blade
LD-5006S-STR
LD-5006C- CVD

Micro Tissue Forceps
Suture Forceps


LDTC-101
Adson
12cm 1x2



LDTTC-102
Adson
12cm



LDTTC-103
Adson, 12cm
Ergonomic Handle



LDTTC-104-15 cm
LDTTC-105-17 cm
Micro Tweezer TC



LDTTC13-400
Micro Tweezer
TC, 17cm



LDTTC13-410
Micro Tweezer
TC, 17cm



LD13-150
DeBakey Atraumatic
Tissue Forceps
Str, 16 cm



LD13-152
DeBakey
Atraumatic
Tissue Forceps
Cvd, 16 cm



LDDDD-102
Corn Suture Pliers
Diamond Coated



LDDDD-100
Corn Suture Pliers
Diamond Coated

DD Coated Point


Medical Grade Titanium Made Instruments



LDTT15-10
Halsey 13cm



LDTT15-17
Mosquito Forceps
12.5 cm, Curved

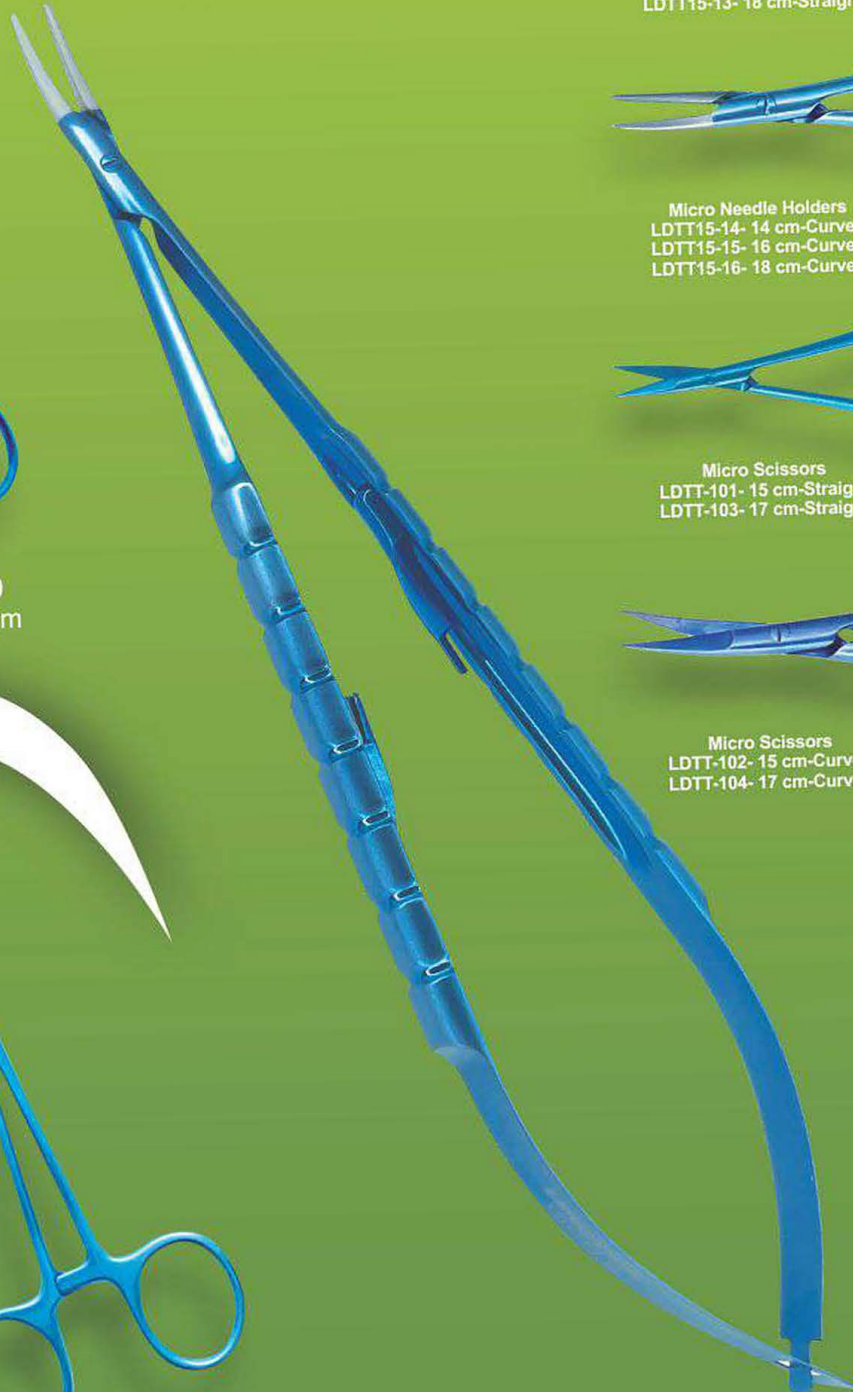

Needle Holders
LDTT15-11- 14 cm-Straight
LDTT15-12- 16 cm-Straight
LDTT15-13- 18 cm-Straight



Micro Needle Holders
LDTT15-14- 14 cm-Curved
LDTT15-15- 16 cm-Curved
LDTT15-16- 18 cm-Curved




Micro Scissors
LDTT-101- 15 cm-Straight
LDTT-103- 17 cm-Straight




Micro Scissors
LDTT-102- 15 cm-Curved
LDTT-104- 17 cm-Curved

LENOX INS Needle holders feature tungsten carbide inserts with a matched diamond-cut pyramidal profile for enhanced longevity, a secure, positive suture needle grasp and reduced risk of cutting the suture material. With a precision designed hinge, the smooth surface reduces the risk of catching and or cutting the suture material during the knot tying process, enhancing clinical efficiencies. For optimal grip and control, LENOX INC offers a variety of design options, including standard, finger and palm controlled closure.

Tungsten Carbide Needle Holders

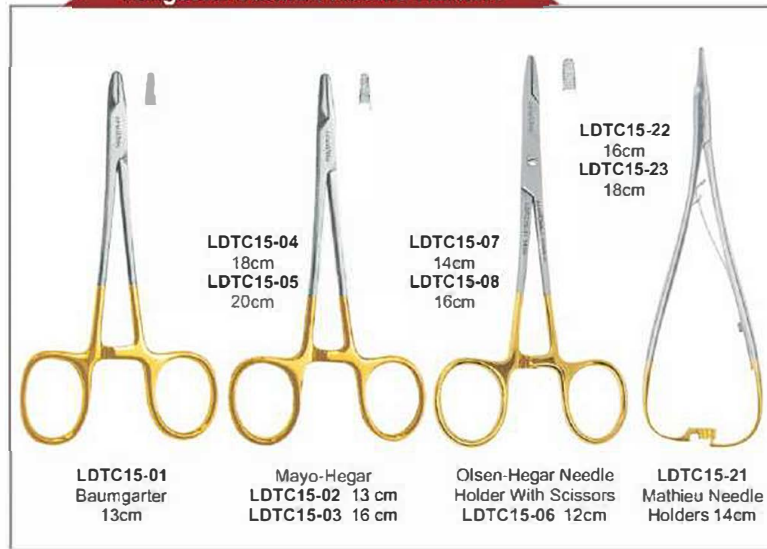

Needle Holders

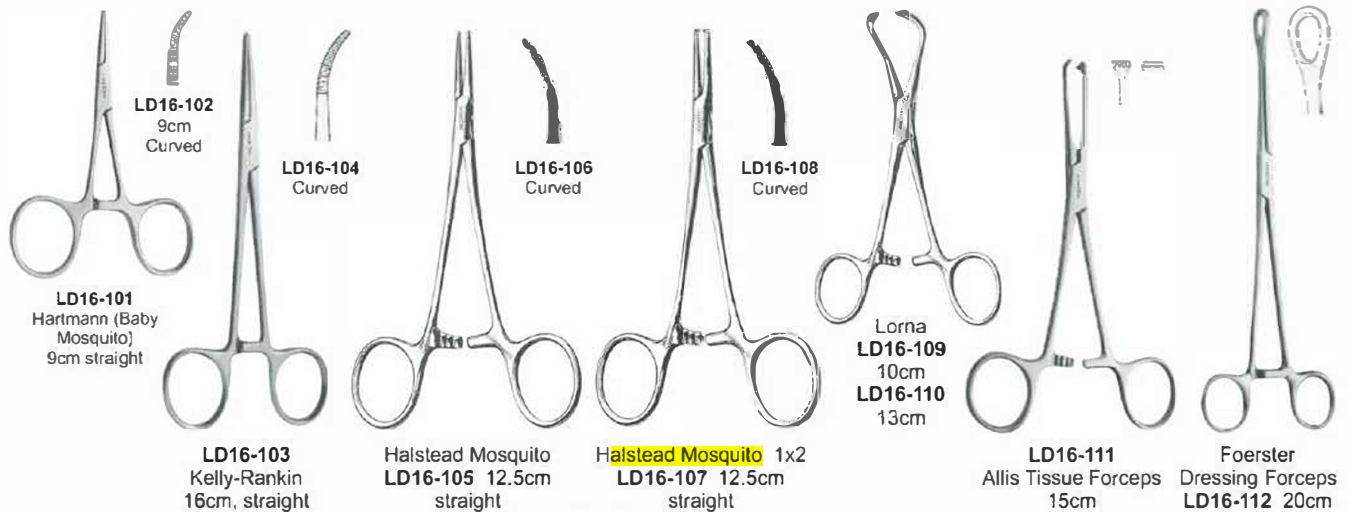

Matrix Holding Instruments





A1-Double Ended



A2-Double Ended




A3-Double Ended



A4-Double Ended



A5-Double Ended

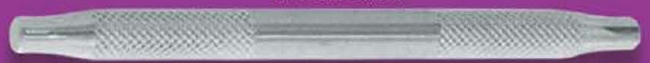

A11-Double Ended



A6-Double Ended



A7-Double Ended



A8-Single Ended



A9-Single Ended



A10-Double Ended



A12-Double Ended



A13-Single Ended

TITANIUM COATED HANDLES (B)

HOLLOW HANDLES (A)



B1-Double Ended



B2-Double Ended



B3-Double Ended



B4-Double Ended



B5-Double Ended



B6-Double Ended




B7-Single Ended




B8-Single Ended



SOLID HANDLES (D)

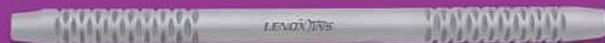


D1-Double Ended

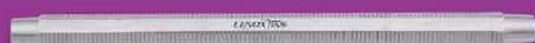

D2-Double Ended



D3-Double Ended

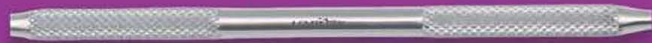

D4-Double Ended



D5-Single Ended



D6-Double Ended

SILICON HANDLES (E)



E1-Double Ended



E2-Double Ended



E3-Double Ended



E4-Double Ended



E5-Double Ended



E6-Double Ended

Hollow Hands (A)



Titanium Coated Hands (B)




Solid Hands (D)




Silicon Hands (E)



Explorers Double Ended

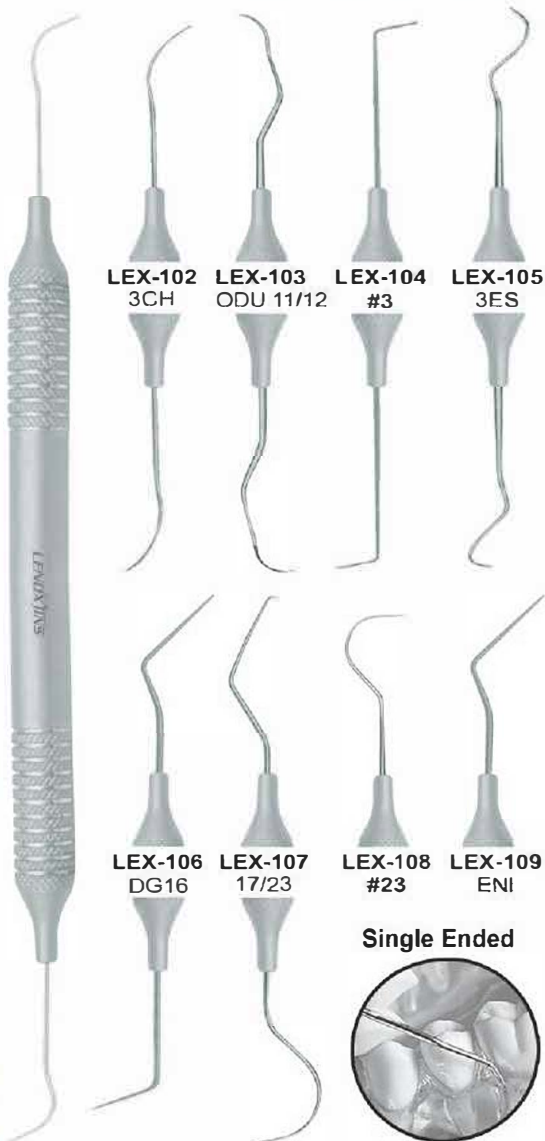

Probes Double Ended

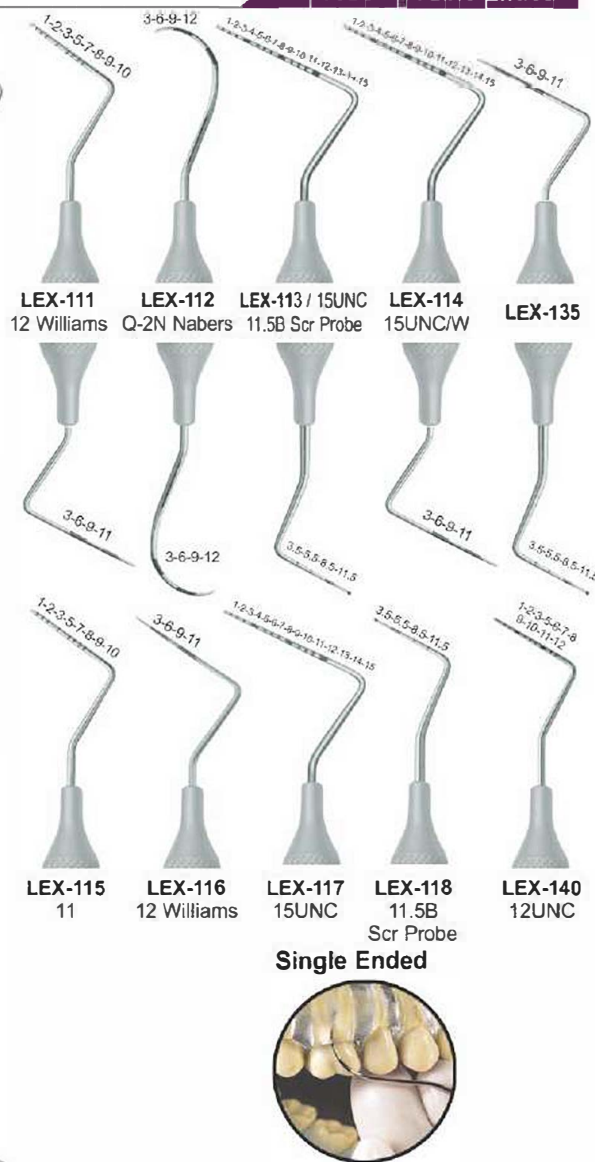


LEX-101
Explorers Double
Ended 2A

LEX-110
Probe Double Ended 12

EXPOS




A double-ended instrument that combines two critical instruments required for every effective diagnostic set-up with an explorer on one end and a probe on the other, an Expos is excellent for dental examinations and charting.

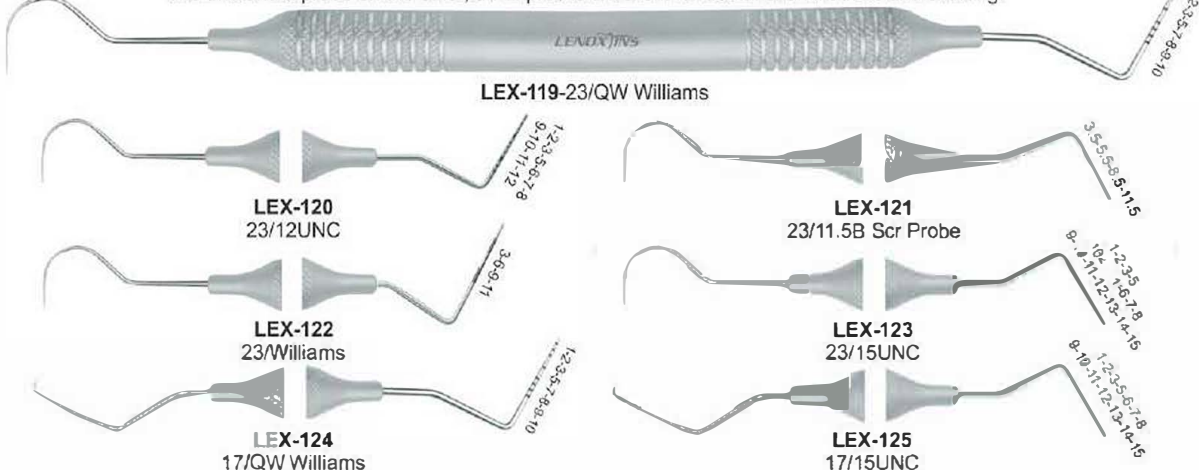

Root Canal Double Ended

Root Canal Single Ended

Hollow Handles (A)

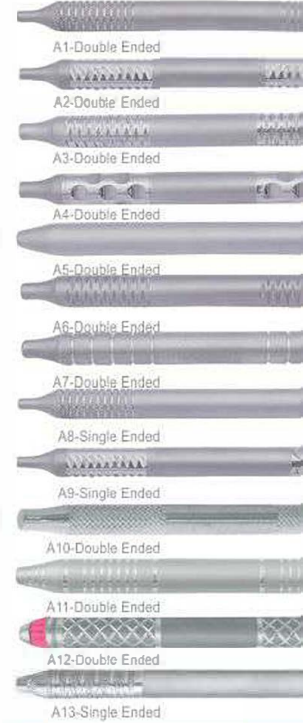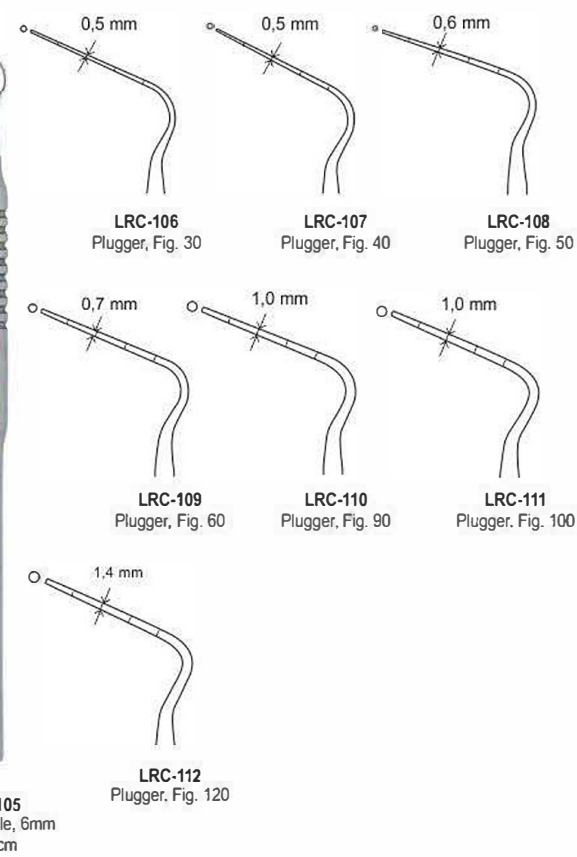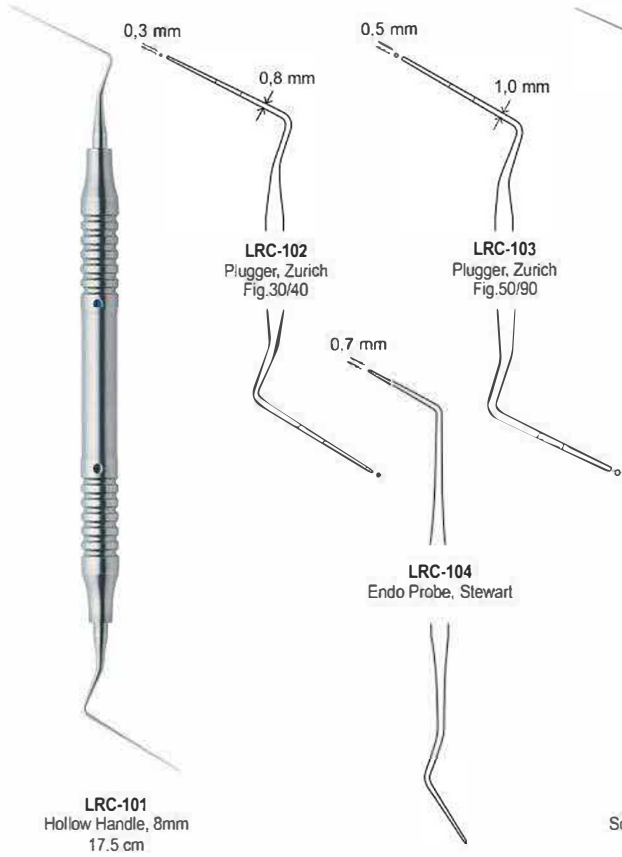


Root Canal Single Ended


Titanium Coated Handles (B)




Solid Handles (D)

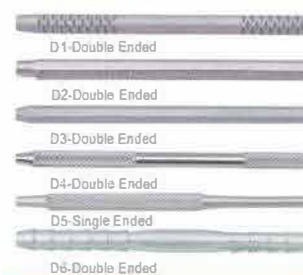

Silicon Handles (E)

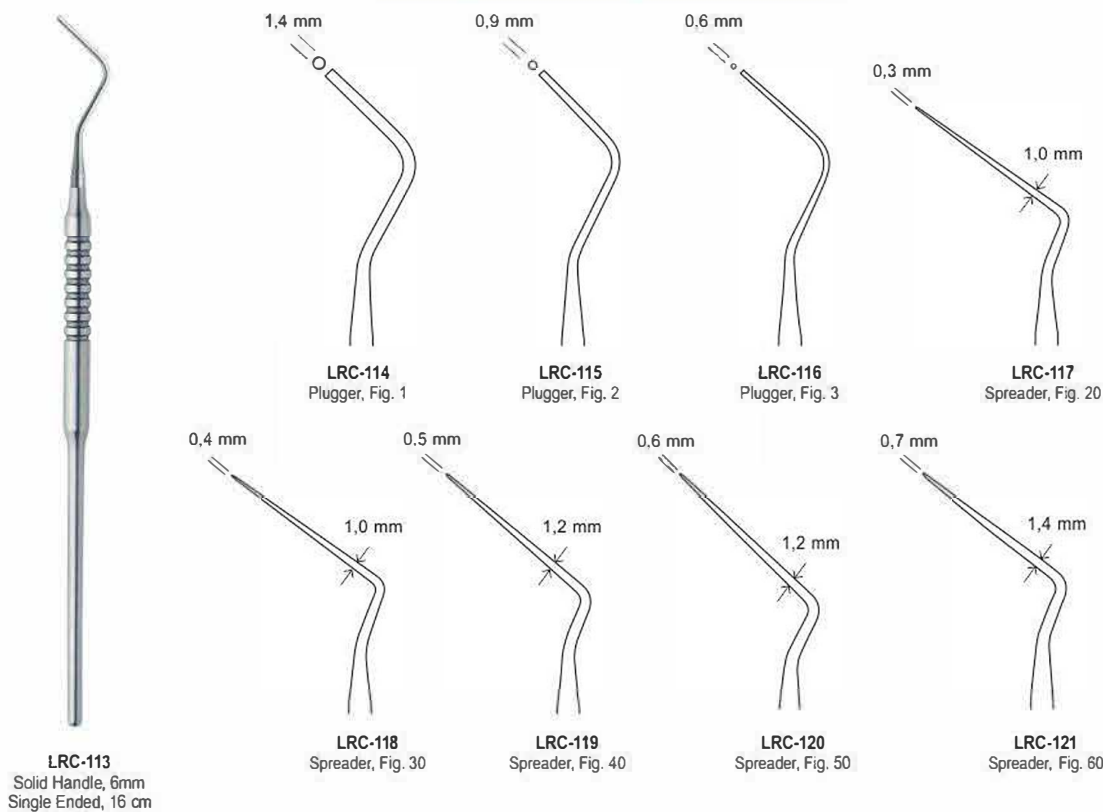

Hollow Hands (A)

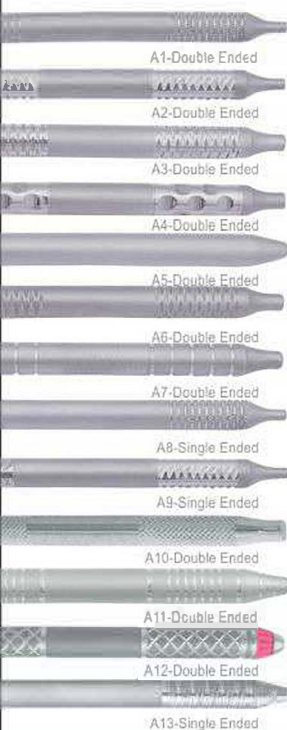

Titanium Coated Hands (B)



Solid Hands (D)




Silicon Hands (E)



Excavators


SPOONS: For Removal of Carious Dentin

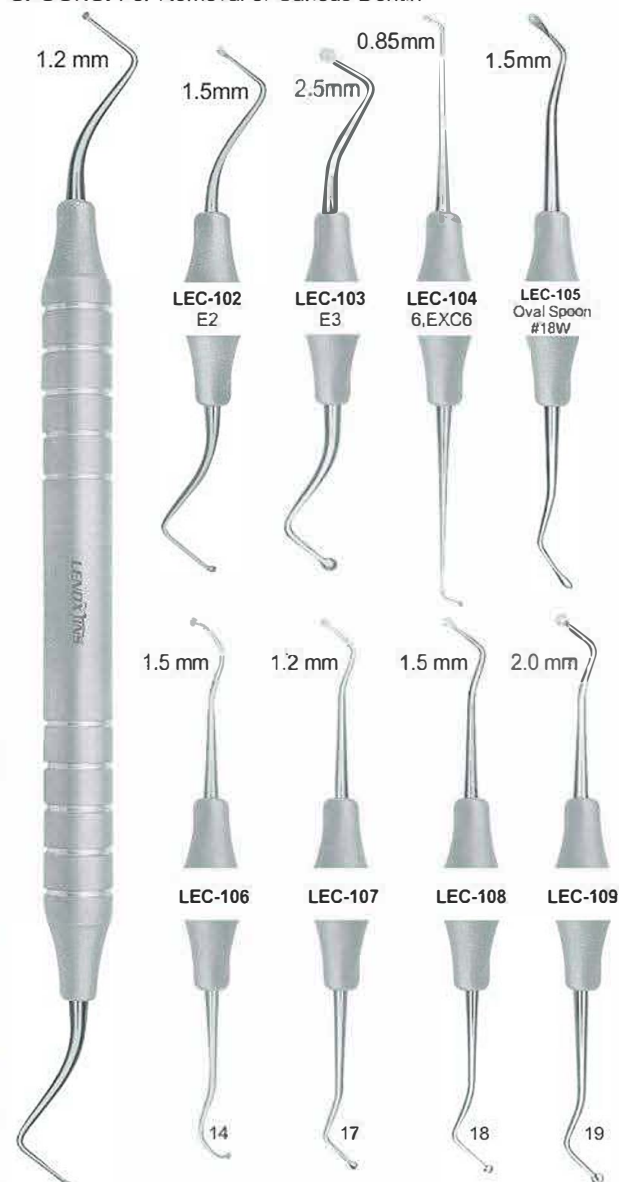

LEC-101-E1

Pluggers / Condensers


Used to Pack Filling Materials into Cavity Preparations

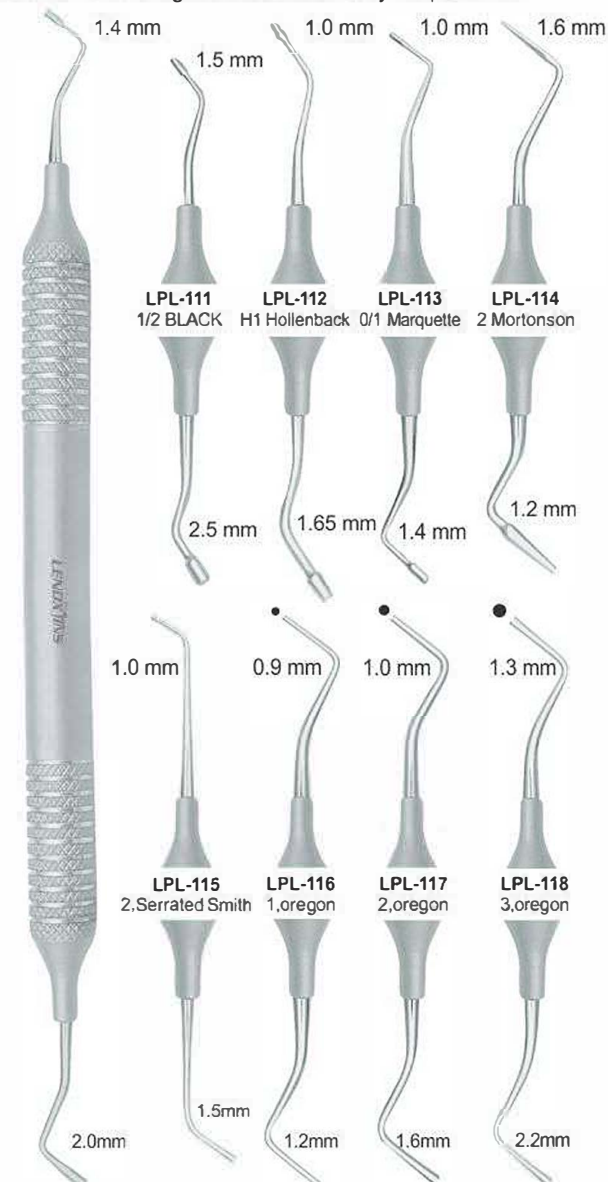

LPL-110-1/2, Andrew


CARVERS

Used to Crave anatomical features and trim excess materials.

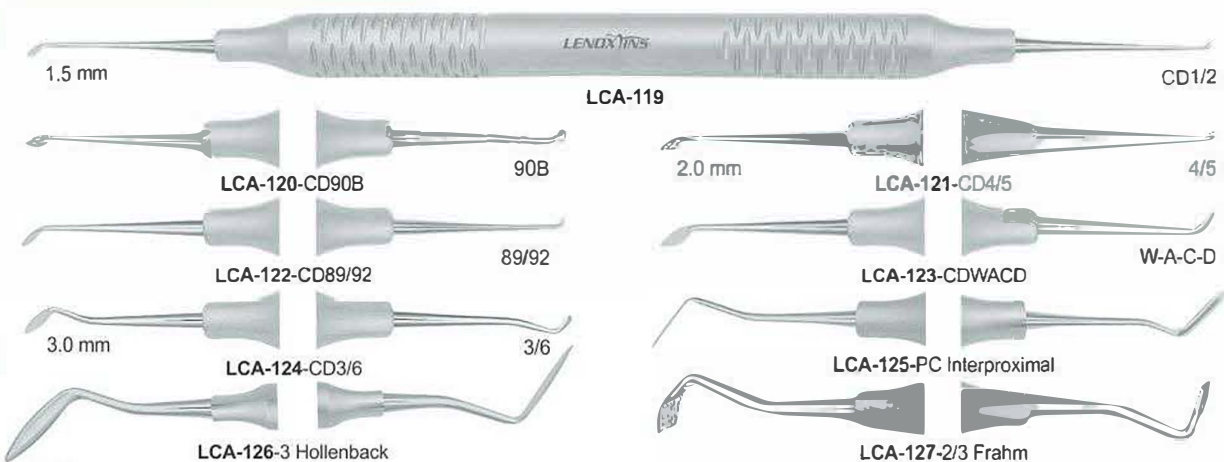

Burnishers

Designed to condense, smooth, carve and polish amalgam.
Gold Titanium Coated Points.

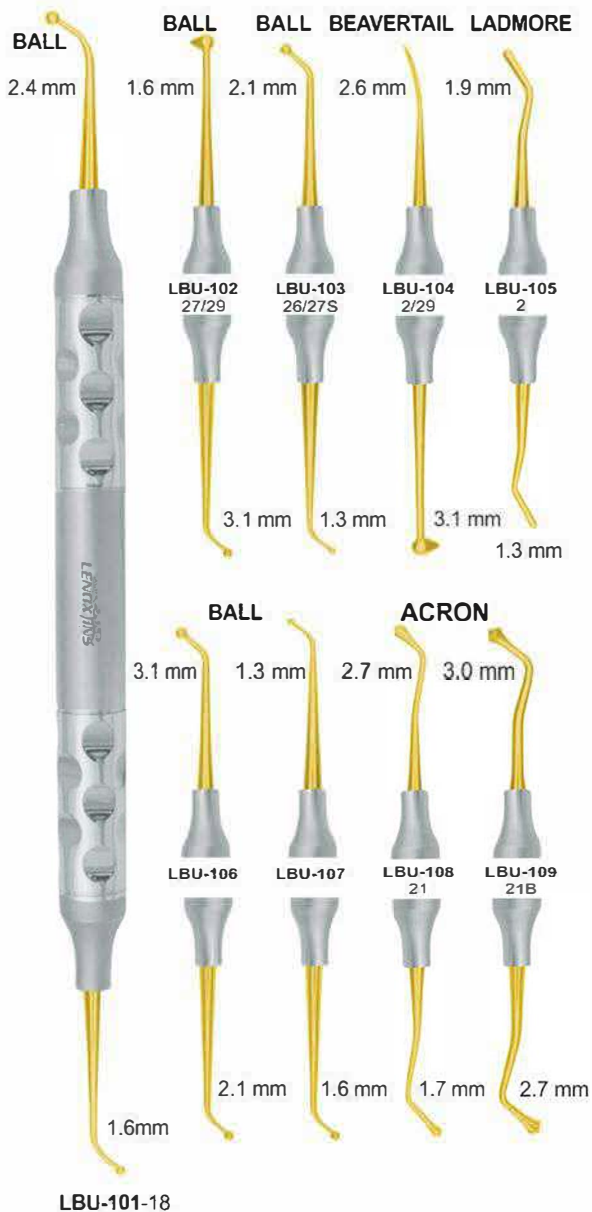

LBU-101-18

Carving Knives

Used for trimming excess filling material, flash and overhangs.



LCK-113 / 20 Esthetic

For anterior teeth. Sharp, thin blade allows access to all surfaces. The offset angle provides universal adaptability.

LCK-114 / 21 Esthetic

For posterior teeth. Thin, sharp offset angle provides access to surfaces where linear finishing strips would not be effective.

LCK-115 / 6 Tanner


For posterior teeth. Sharp offset angle provides access to many surfaces.

LCK-116 / 12 Novatech Cement Remover

Combines a sickle-shaped scaler with a flat blade for removal of excess resin, cement, or porcelain flash.
The narrow chisel removes excess interproximal material with a push stroke.

PK-Thomas

Used for waxing procedures and techniques.

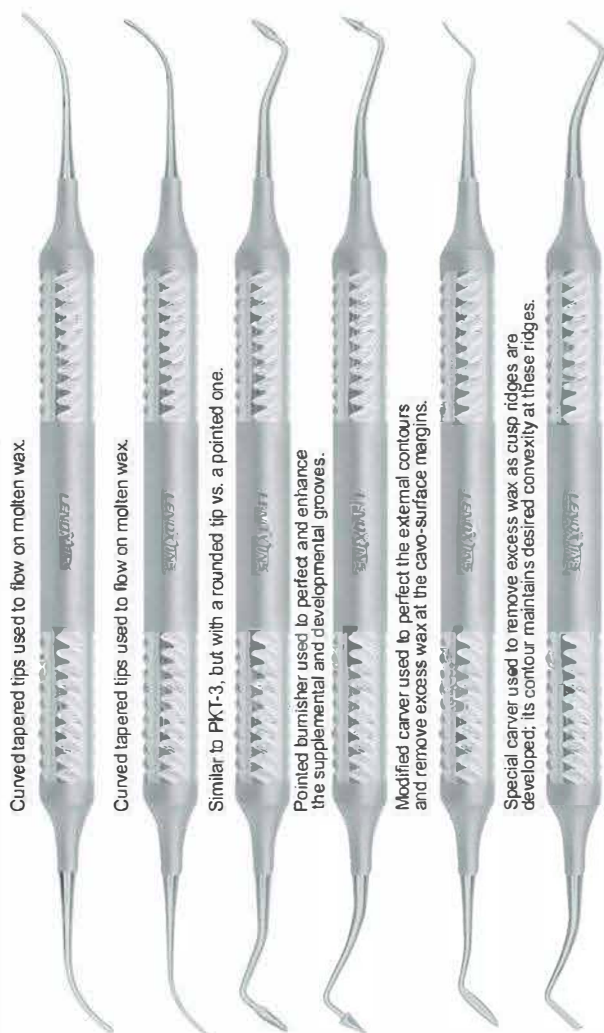

LPKT-101 LPKT-102 LPKT-103 LPKT-103R LPKT-104 LPKT-105

Hollow Hand's (A)



Titanium Coated Hands (B)



Solid Hands (D)



Silicon Hands (E)



Hollow Hands (A)

A1-Double Ended

A2-Double Ended

A3-Double Ended

A4-Double Ended

A5-Double Ended

A6-Double Ended

A7-Double Ended

A8-Single Ended

A9-Single Ended

A10-Double Ended

A11-Double Ended

A12-Double Ended

A13-Single Ended

Titanium Coated Hands (B)

B1-Double Ended

B2-Double Ended

B3-Single Ended

B4-Double Ended

B5-Double Ended

B6-Double Ended

B7-Single Ended

B8-Single Ended

Solid Hands (D)

D1-Double Ended

D2-Double Ended

D3-Double Ended

D4-Double Ended

D5-Single Ended

D6-Double Ended

Silicon Hands (E)

Cord Packing Instruments

For atraumatic and accurate cord placement.

LCP-116 / BN1

Thin blade and rounded contour facilitates use in both thick and thin tissues without catching or dropping cord. Bilateral notch allows placement in limited access areas.

LCP-117 / 1 Yardley

Round non-serrated working end.

Serrated

LCP-118
CSI-1

Non-Serrated

LCP-119
7 Guyer

Serrated

LCP-120
113

Non-Serrated

Gingival Retractors

Protects tissue during cavity preparations and placement of composite restorations.

LGR-120 / 1 Meinershagen / For mandibular incisors

LGR-121 / 2 Meinershagen

For maxillary and mandibular premolars and canines.

LGR-122 / 3 Meinershagen


For maxillary central incisors and wide canines

LGR-123 / 4 Meinershagen
For all molars.LGR-124
10 Goldman-Fox


Composite Filling Instruments

Also Available in Titanium Coated Colored Points.


Vertical and horizontal combination.

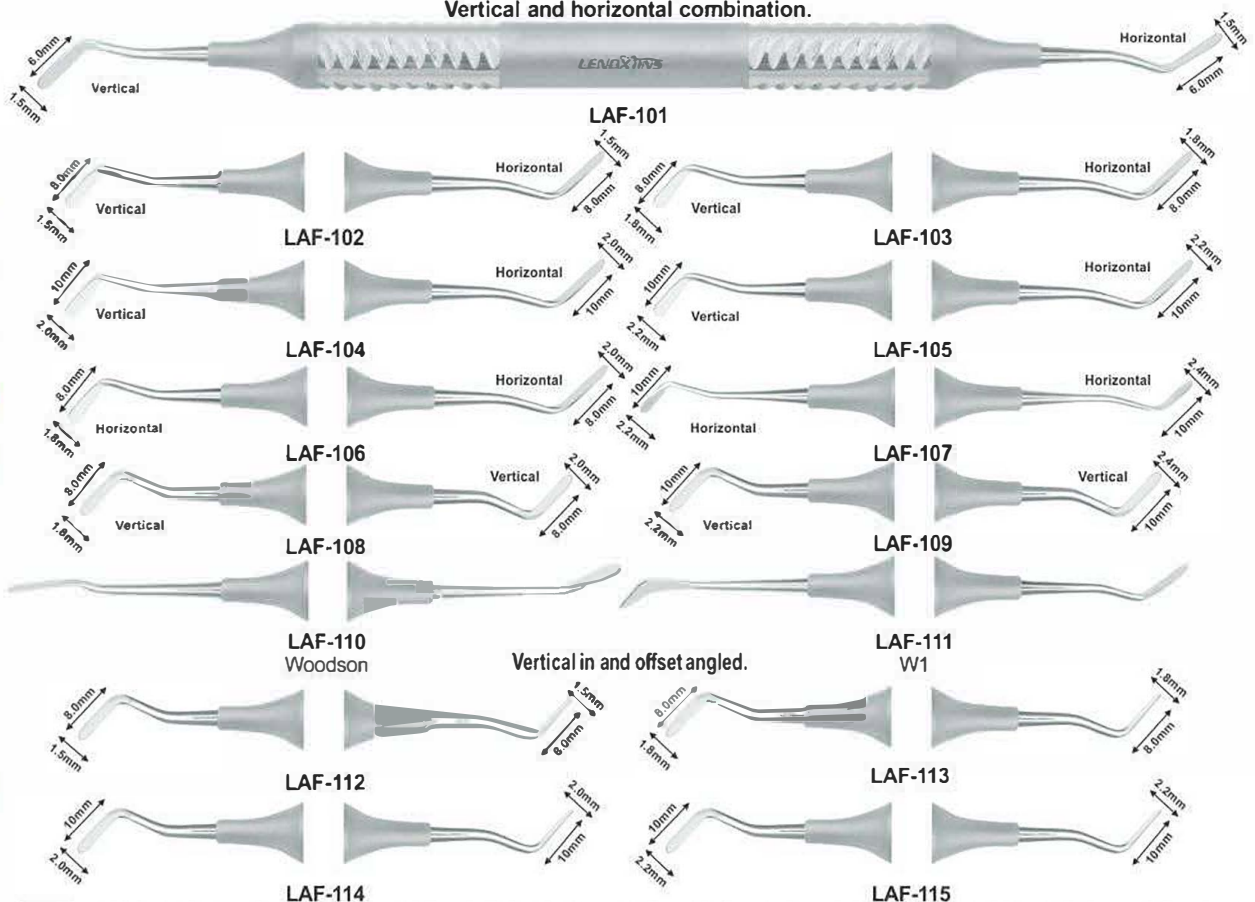


Composite Filling Instruments

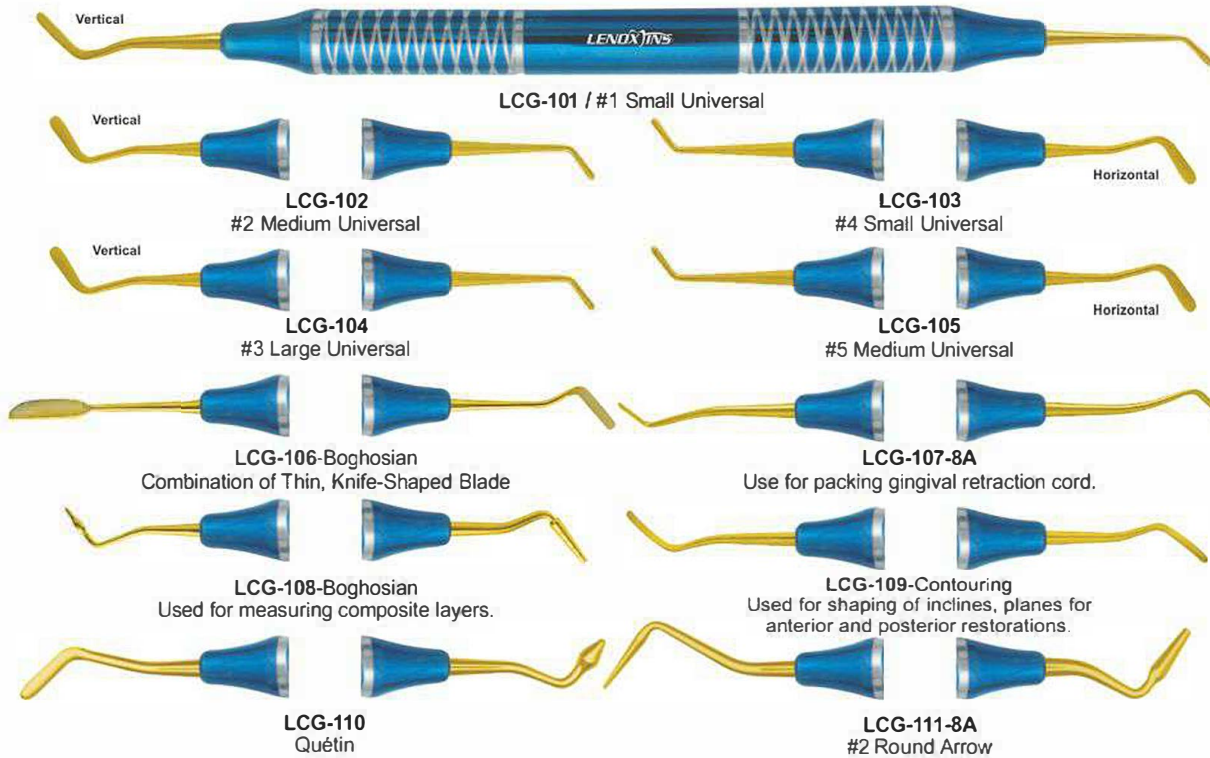

Composite Filling Instruments Set

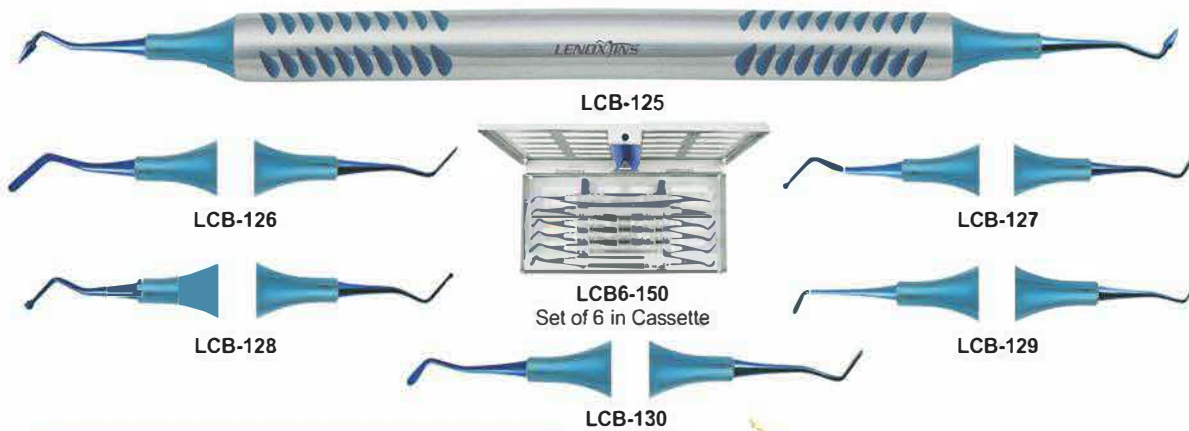

Amalgam Carrier

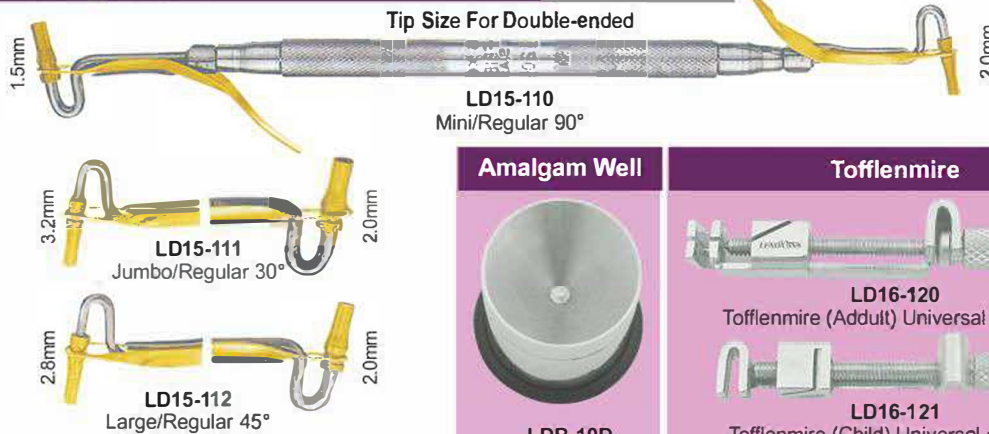

Amalgam Well



Tofflemire



Hollow Hands (A)



Titanium Coated Hands (B)



Solid Hands (D)



Silicon Hands (E)



Hollow Hands (A)

A1-Double Ended

A2-Double Ended

A3-Double Ended

A4-Double Ended

A5-Double Ended

A6-Double Ended

A7-Double Ended

A8-Single Ended

A9-Single Ended

A10-Double Ended

A11-Double Ended

A12-Double Ended

A13-Single Ended

Titanium Coated Hands (B)

B1-Double Ended

B2-Double Ended

B3-Single Ended

B4-Double Ended

B5-Double Ended

B6-Double Ended

B7-Single Ended

B8-Single Ended

Solid Hands (D)

D1-Double Ended

D2-Double Ended

D3-Double Ended

D4-Double Ended

D5-Single Ended

D6-Double Ended

Silicon Hands (E)

E-6

E-5

E-4

E-3

E-2

E-1

Sickle Scalers

Designed for removal of supragingival calculus.

All Scalers & Cutters are Available in

J2-420 ECO

Best Quality Medical Stainless Steel.
High Chromium Content
Good Strength
Easy To Sharpe
Economical Cost
Available in Different Sizes.

USA 440c Premium

440c High-end Stainless Steel
High Chromium Content
High Carbon Martensitic
High Strength
Better Edge Retention
Titanium Coated Points



LSC-101

Contra angle design offers greater access to distal surfaces



LSC-102 / H6/H7

Super-thin, contra angle design is for interproximal areas.
Shank angulation provides accessibility to all anterior and premolar surfaces



LSC-104 / H5/33 - Hygienist-Jacquette

Super-thin sickle scaler combined with a straight sickle blade.
For scaling interproximal and cervical areas of anterior teeth.

LSC-106
U15/33

Towner-Jacquette




LSC-108

For posterior supra-gingival
interproximal scaling



LSC-110/204SD

Contra angle design offers greater
access to distal surfaces

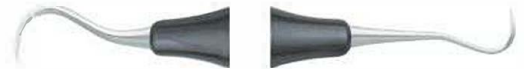LSC-112
Nabrasca 5/6LSC-114
Nevi-1 / H5
AnteriorLSC-103
H5/U15
Hygienist-
TownerLSC-105
U15/30
Towner-JacquetteLSC-107
204

LSC-109 - 0/00 Morse

Miniature blade for scaling crowded or overlapping
anterior teeth excellent for stain removal

LSC-111
204SLSC-113
NV1/H5LSC-115
Nevi-2 / Posterior
A super thin posterior curved sickle scaler

Columbia Curettes

LCC-106
13/14 Columbia
(Regular)

Universal Application

LCC-107
13/14 Columbia
RigidLCC-108
2R/2L Columbia, Anterior applicationLCC-109
4R/4L Columbia, Posterior application


Langer Currettes



Barnhart




Periodontal Files

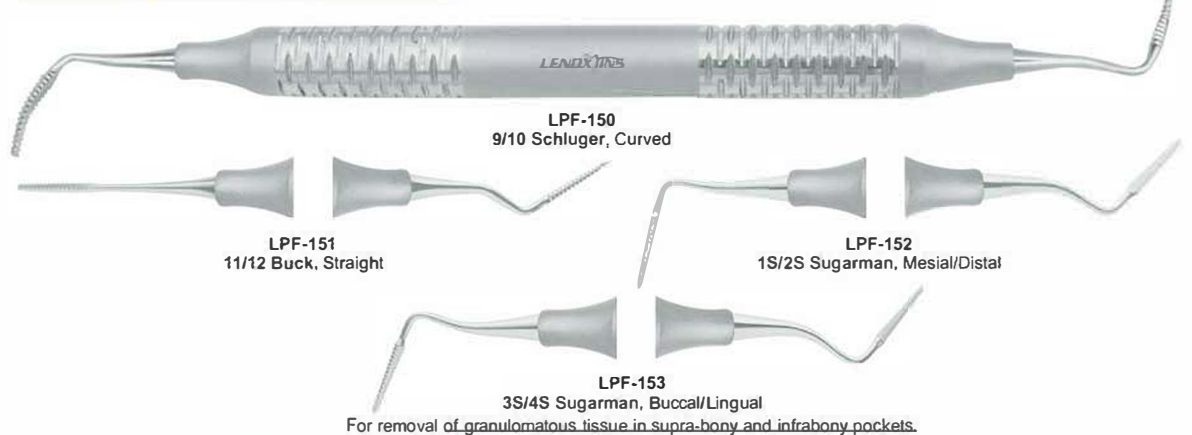

Goldman-Fox



Hollow Handls (A)

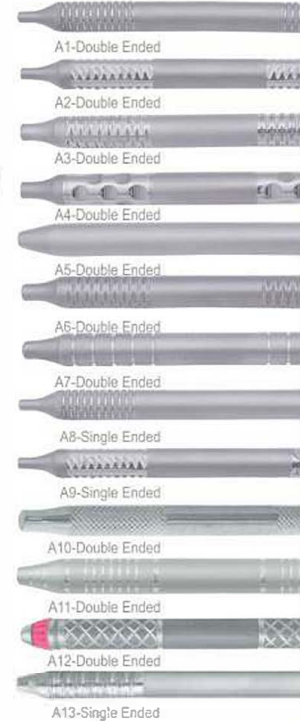

Titanium Coated Handls (B)



Solid Handls (D)



Silicon Handls (E)




- Area-specific to allow for deep scaling, root planing and periodontal debridement
- Offset blade provides a perfect working angulation for the tooth surface
- Offset blade angulation means only the lower cutting edge of each blade is used

Gracey Standard Curettes




Hollow Hands (A)

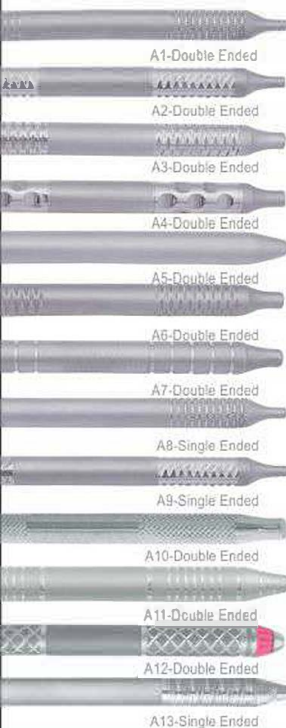

Titanium Coated Hands (B)



Solid Hands (D)

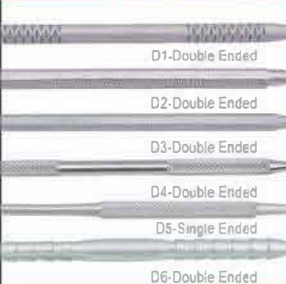


Silicon Hands (E)



		Slight contra-angle for anterior incisors and canines		
LGC-110 1/2 (Reg-shown)	LGC-110-1 1/2 (Rigid)		LGC-110-2 1/2 After Five	LGC-110-3 1/2 Mini Five
		Short contra-angle for anterior incisors and canines		
LGC-111 3/4 (Reg-shown)	LGC-111-1 3/4 (Rigid)		LGC-111-2 3/4 After Five	LGC-111-3 3/4 Mini Five
		Medium contra-angle for anteriors and premolars		
LGC-112 5/6 (Reg-shown)	LGC-112-1 5/6 (Rigid)		LGC-112-2 5/6 After Five	LGC-112-3 5/6 Mini Five
		Medium contra-angle for premolars and molars, facial and lingual surfaces		
LGC-113 7/8 (Reg-shown)	LGC-113-1 7/8 (Rigid)		LGC-113-2 7/8 After Five	LGC-113-3 7/8 Mini Five
	Long contra-angle for molar buccal and lingual surfaces and a more pronounced angle than 7/8 for less accessible root surfaces.			
LGC-114 9/10 (Reg-shown)			LGC-114-1 9/10 (Rigid)	
		Angulated to reach mesial surfaces of posterior teeth		
LGC-115 11/12 (Reg-shown)	LGC-115-1 11/12 (Rigid)		LGC-115-2 11/12 After Five	LGC-115-3 11/12 Mini Five
		Angulated to reach distal surfaces of posterior teeth		
LGC-116 13/14 (Reg-shown)	LGC-116-1 13/14 (Rigid)		LGC-116-2 13/14 After Five	LGC-116-3 13/14 Mini Five
		Same shank bends as 13/14, but blade positioned to reach mesial posterior surfaces		
LGC-117 15/16 (Reg-shown)	LGC-117-1 15/16 (Rigid)		LGC-117-2 15/16 After Five	LGC-117-3 15/16 Mini Five
	Accentuated angles improve access to distal posterior surfaces. The 17/18 also has smaller blade and slightly longer terminal shank.			
LGC-118 17/18 (Reg-shown)			LGC-118-1 17/18 (Rigid)	

Periosteals

For reflecting and retracting the mucoperiosteum after incisions of the gingival tissue.

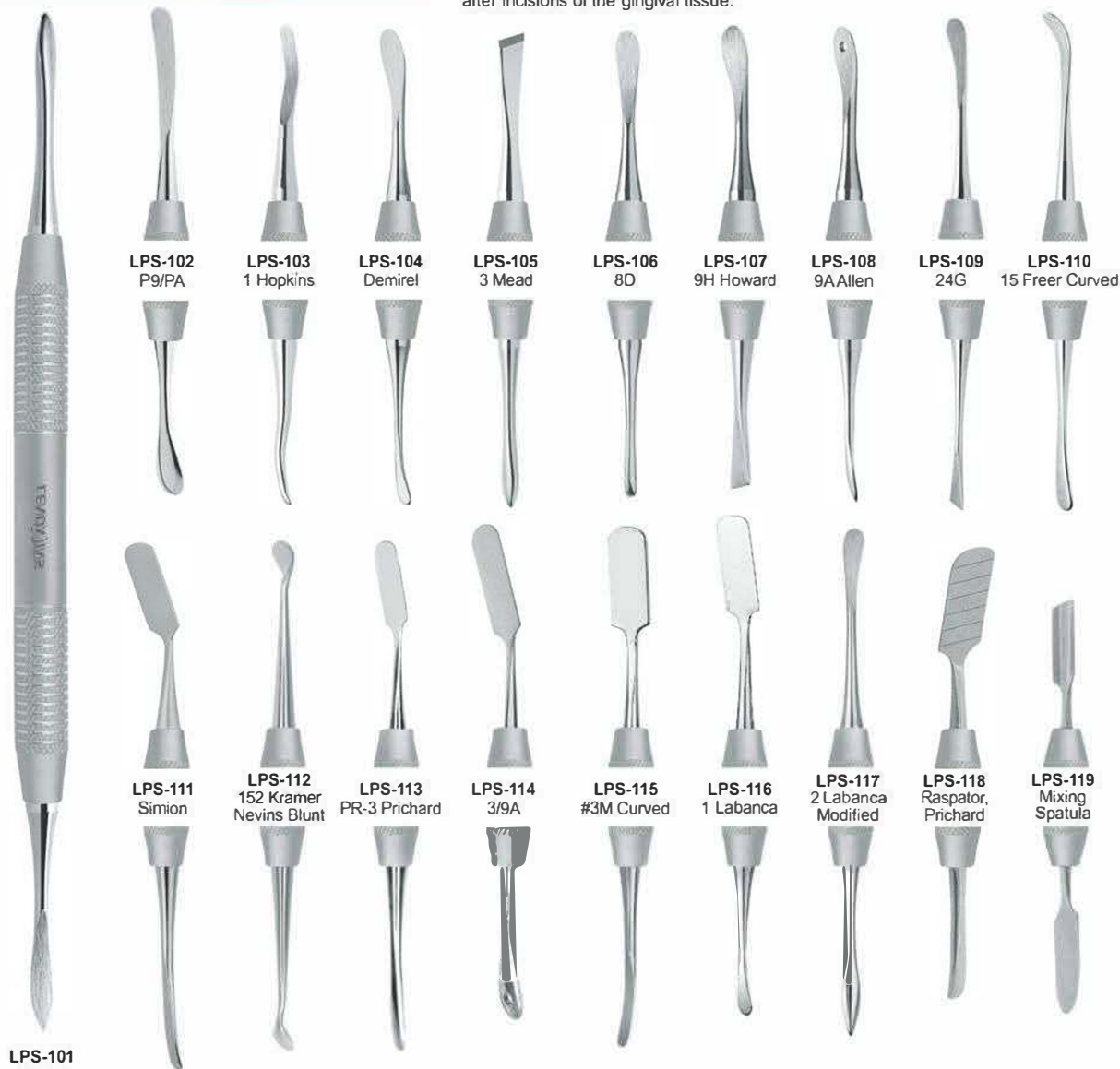

Periosteal Elevators



Hollow Handls (A)



Titanium Coated Handls (B)



Solid Handls (D)



Silicon Handls (E)



Retractors & Mouth Gags



LD17-101-Infant, Small
LD17-102-Children, Medium
LD17-103-Adults, Large



Weitlaner
LD17-104-10.5cm
LD17-105-13.5cm



LD17-106
Stainless Steel
Double Ended Cheek
Retractors



LD17-107
Stainless Steel
Single Ended Cheek
Retractors



Cheek Retractors,
LD17-108-14cm
LD17-109-16 cm



Minnesota
LD17-110-Small
LD17-111-Large



LD17-112
Retractors
17 cm
Hooks 22 mm



LD17-113
Retractors, Cheek
retractors combination
21 cm, 2



LD17-114
140 mm
LD17-115
150 mm



LD17-116
Bishop



LD17-117
Austin Fig.1



LD17-118
Austin Fig.2

Surgical Retractor
LD17-119-7x25mm
LD17-120-10x42mm
LD17-121-12x55mm

LD17-122
Implant
Retractor
23cm

LENOX INS Adjustable Scalpel Handles



Locks At 90°, 60°, 30°
Straight, 120°, 150°, 180°
Angles

LDAS-1
6 Way Adjustable Angled
Scalpel Handle

Scalpel Handles



LD18-101 No. 5E LD18-102 No. 5A LD18-103 3K Mini LD18-104 Fig.3 LD18-105 Fig.4 LD18-106 Fig.3G

New Scalpel Handles




LD18-110 Double Blade holder LD18-120 Push Scalpel Handle

Cement Spatulas & Wax Knives



LD20-101 Gritman LD20-102 3 LD20-103 1 Medium Cement LD20-104 Small Wooden Handle LD20-105 Large Wooden Handle LD20-106 Wax Knives, 20 cm curved

Rubber Dam Forceps

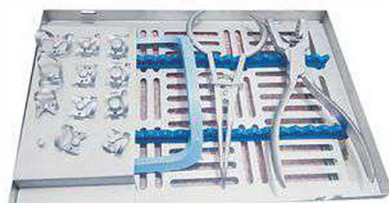

LDS21-100
01-Anisworth
02-Brewer
03-Clamps Stand
04-Rubber Dam Frame
05-Rubber Dam Clamp
Fig.212,211,8,8A
14,0,00,2,2A,7



LDS21-101
01-Anisworth
02-Palmer
03-Clamps Stand
04-Rubber Dam Frame
05-Rubber Dam Guide
06-Rubber Dam Clamp
Fig.212,211,8,8A
14,0,00,2,2A,7



LD21-102
Anisworth



LD21-103
Brewer



LD21-104
Ivory light weight



LD21-105
Palmer



LD21-106
Seaparating Elastic Plier



LD21-170
Dental Ring Applicator

Rubber Dam Clamps



37



LD22-101
Dental
crown remover
forceps



LD22-102
Crown Spreaders



LD22-103
WYNNAN
Crown Gripper, CARB-N-SERT



LD22-104
Key Crown Remover
Set of 3 pcs



LD22-105
Universal
Crown Remover
Hollow Handle



LD22-106
Automatic
Crown Remover
with 2 Tips



LD22-100
Crown Remover Set
Full Automatic



LD23-100
Boley Gauge



LD22-107
Napkin Holder
Silicon

LD22-108
Napkin Holder

Bone Gauges and Measuring Instruments

Bracket Positioning Boons Gauge

LD23-101

- It is The Most Frequently Used Bracket Positioning Aid.
- It is Made Of Stainless Steel.
- It Accurately Measures Height Of Bracket Placement From The Incisal Edge At 3.5, 4, 4.5 And 5mm.
- Flat Surface Rests On Incisal/occlusal Surface Of Tooth. Simplify Seating The Bracket On The Teeth With Exactly The Right Distance Between Bracket Dot And Incisal Edge With Pencil Lid Or Meial Tips.



LD23-102



LD23-103



LD23-104
Measuring Loop
Straight




LD23-105
Measuring Loop
Curved

LD23-106
Lwanson Crown
Gauge

Bone Pluggers



Sinus Lift Instruments

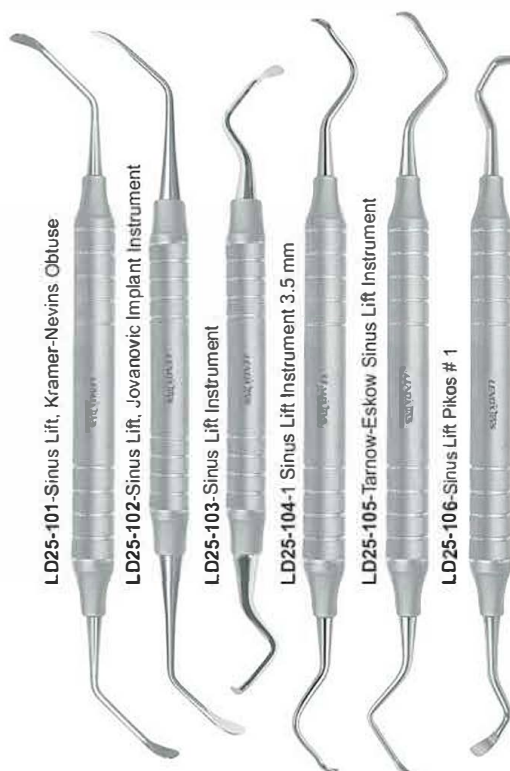


Bone Chisels



Tissue Punch Straight & Curved



Convexes Straight & Curved

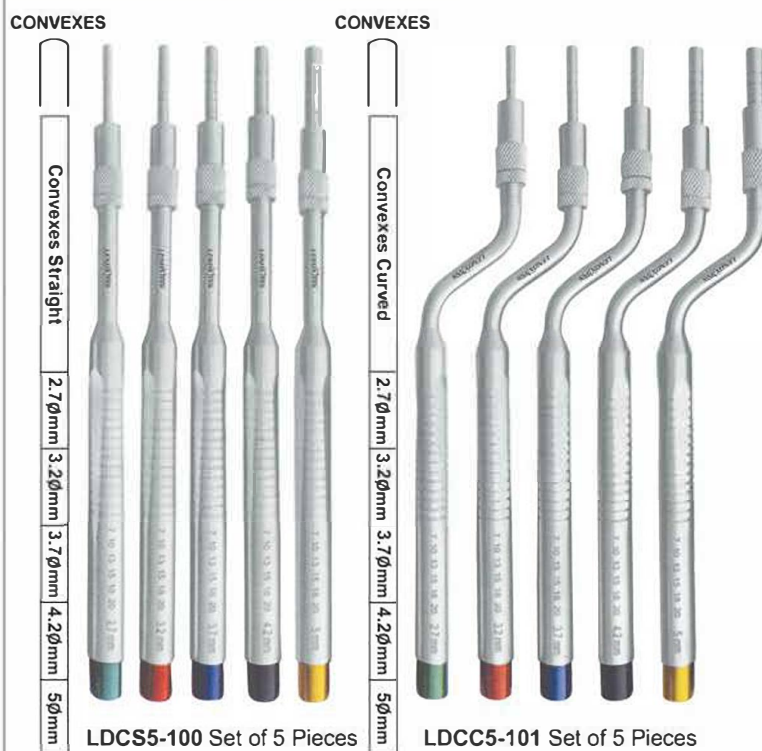



LD29-101
Bone Ridge Mapping
Forceps



LD29-102
Measuring Caliper



LD29-103
Splinter Forceps
Straight



LD29-104
Splinter
Forceps
Curved



LD29-105
Bone Mill Plier



LD29-106
Molt 18.5cm
Auto Clavable
Plastic Tips.

Bone Scraper



LD30-100
Curved



LD30-101
Straight

Star Torque Wrench



LD31-100
Star Torque Wrench
6.35mm Universal



LD29-108
Bone Grafting Syringe
17 cm, Straight / Bend
I.D 3.5mm, O.D 4.5



LD29-109
Bone Grafting Syringe
18 cm, Straight / Bend
I.D 4.5mm, O.D 5.5



LD29-110
Bone Grafting Syringe
18 cm, Straight / Bend
I.D 6.0mm, O.D 7.0



LD29-107
Bone Grafting Syringe
17 cm, Straight / Bend
I.D 2.5mm, O.D 3.5



LD29-111
Suction
Tub



LD29-112
Fig. 6FR
Frazier
15 cm-Ø2mm



LD29-113
Ø3mm
Length 180mm



LD29-114
Ø4mm
Length 180mm




LD29-115
Ø5mm
Length 180mm



LD29-130
VDO GAUGE
20-100mm



LD29-132
20-130mm

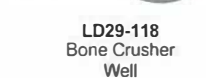

LD29-118
Bone Crusher
Well

Bone Expander Kit



LDBE-100
Bone Expander Kit
Expander: 2.6mm, 3.0mm
3.4mm, 3.8mm, 4.3mm
Ratchet Driver
Parallel Pin
Pilot Drill
Drill 1.4mm
Drill Driver 4mm
Hand Driver 4mm
Ratchet 4.0mm

Periotomes



LD32-100 / P15 Titanium TG CVD S/E



LD32-101 / P14 Titanium um STR S/E



LD32-102 / P13 Titanium Angle S/E



LD32-103 / P11 Titanium AV Regular



LD32-104 / P12 Titanium AV Regular

Micro Serrated





LD32-120
Mini Drill Set of 6 Pcs.
in Stainless Steel Stand,



LD32-121
Standard Drill Set of 6 Pcs.
in Stainless Steel Stand,



LD32-122
Trephine Bur Set of 6 Pcs.
in Stainless Steel Stand,



LD32-123
Bone Compression
Prosthodontic Dental Kit
for Implant



LD32-124
Tissue Punch Kit
in Stainless Steel Stand



LD32-125
Trimmer Kit in
Stainless Steel Cassette,




LD32-126
Mini Expander Kit
in Stainless Steel Cassette



LD32-127
Bur Holder
for 12 Pcs

LD32-128
Bur Holder
for 06 Pcs

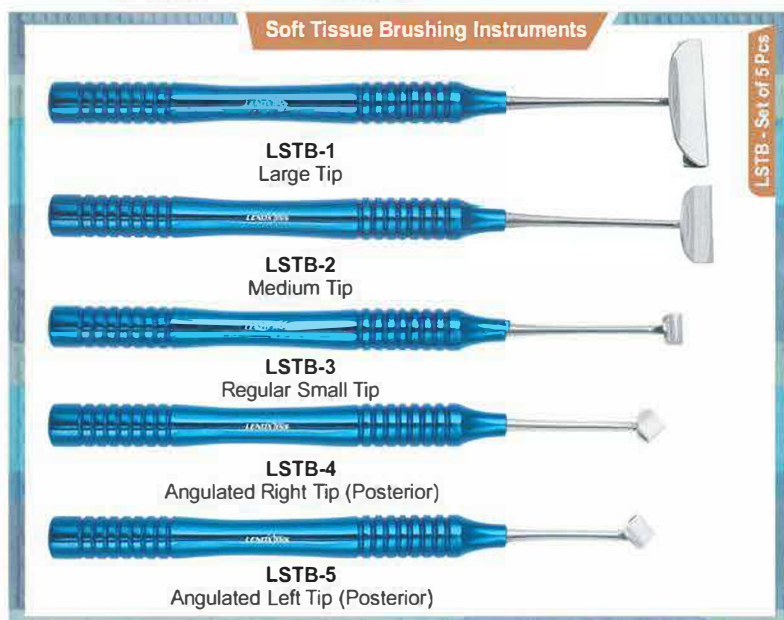

LD32-129
Conical Universal Kit
(Kit of 50 Pcs)
in Stainless Steel Cassette



Distal End Cutters



Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)

LD33-101
Distal End
Cutter Universal cut
long handle 14cm



LD33-102
Distal End Cutter
Universal cut
long Neck 12cm
Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)



LD33-103
Distal End Cutter
Universal cut
long handles 14cm
Max: .021"(.54mm)
x .025"(.64mm)



LD33-104
Flush Cutter with
Hold 12cm
Max: .021"(.54mm)
x .028"(.71mm)
Min: .012"(.31mm)



LD33-105
Distal End Cutter
Universal cut
long handle 14.5cm
Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)



LD33-106
Distal End Cutter
Universal cut
Box Joint 12cm
Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)



LD33-107
Distal End Cutter
Universal cut
long neck 12cm
Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)



LD33-108
Distal End Cutter
Universal cut heavy
pattern.Box Joint long
Handles 14cm
Max: .021"(.54mm)
x .025"(.64mm)
Min: .012"(.31mm)

Ligature & Pin Cutters



Max: .012"(.31mm)
Soft Wires

LD33-109
Ligature & Pin Cutters
straight 12.5cm



LD33-110
Mini Ligature & Pin
Cutters 15 degree
angled, Slim Beak
12.5cm
Max: .012"(.31mm)
Soft Wires



LD33-111
Ligature & Pin
Cutters wide beak
15 degree angled
12.5cm
Max: .012"(.31mm)
Soft Wires



LD33-112
Ligature & Pin
Cutters straight
12cm
Max: .012"(.31mm)
Soft Wires



LD33-113
Micro Ligature & Pin
Cutters 15 degree
angled 11.5cm
Max: .012"(.31mm)
Soft Wires



LD33-114
Ligature Cutters
45 degree angled
wide Handles
very Fine Tips.
11.5cm
Max: .012"(.31mm)
Min: .008"(.20mm)
Soft Wires



LD33-115
Ligature & Pin
Cutters, Box
Joint 1/3 Handles
Gold Platted 12.5cm
Max: .012"(.31mm)
Soft Wires



LD33-116
Micro Mini Ligature
Cutters 15 degree
angle 11cm
Max: .012"(.31mm)
Soft Wires

Pliers



LD33-117
Bracket Removing
Plier 12.5cm angled



LD33-118
Bracket Removing
Plier 12.5cm
straight



LD33-119
Band Removing
Pliers long neck
12.5cm



LD33-120
Slim Weingart
Plier 13cm



LD33-121
Weingart
Plier 12.5cm



LD33-122
Weingart
Plier 14cm



LD33-123
Nance Pliers
12.5cm



LD33-124
Four Step
Loop Forming
Pliers 13cm



LD33-125
Three Jaw Pliers
12.5cm

Max: .028"(.71mm)
Wires

LD33-126
Adre Pliers
13cm

Max: .030"(.76mm)
Wires

LD33-127
Bird Beak
Pliers 12cm

Max: .030"(.76mm)
Wires

LD33-128
Lab Bird Beak
Pliers with serrated
tips 14cm

Max: .036"(.92mm)
Hard Wires

LD33-129
Light Wire Bending
Pliers long 12.5cm

Max: .020"(.51mm)
Wires

LD33-130
Slim Bird Beak
Pliers Tapered
Beak 12.5cm

Max: .022"x0.25"
(.56mm x .64mm)
Wires

LD33-131
Rectangular Arch
Forming pliers
12cm

Max: .022"x0.25"
(.51mm x .64mm)
Wires

LD33-132
Slim Arch
Forming Pliers
12cm

Max: .022"x0.25"
(.51mm x .64mm)
Wires

LD33-133
Tweed Arch
Forming Pliers
12cm

Max: .022"x0.25"
(.51mm x .64mm)
Wires



LD33-134-A
Hollow chop
Arch Forming
Pliers 12cm

LD33-134-B
Hollow chop
Arch Forming
Pliers with 3
Grooves 12cm

LD33-135
Tweed Plier
12.5cm

LD33-136
Adhesive Removing
Pliers

LD33-137
How Utility
Pliers angled
12.5cm

LD33-138
How Utility
Pliers straight
12.5cm

LD33-139
Tweed Pliers
12.5cm

Max: .022"x0.28"
(.56mm x .71mm)
Wires

LD33-140
Loop Forming
Pliers with cutters.
Tweed Type

Max: .025"
Rectangular Hard wires

LD33-141
Step Forming Pliers
A-0.25 / B-0.50
C-0.75 / D-1mm

Max: .022"
(.56mm) Wires



LD33-142
Ligature Forming
Pliers

LD33-143
Niti Cinch Back
Pliers, Slim Beak
Max: .021"x0.25"
(.54mm x .64mm)
Wires

LD33-144
Hook Cimping
Pliers

LD33-145
Loop Forming
Pliers, long Beak
Max: .028" (.70mm)
Hard Wires

LD33-146
Optical Pliers

LD33-147
Young Loop
Bending Pliers
Max: .028" (.70mm)
Hard Wires

LD33-148
Crown Crimping
Pliers

LD33-149
Johnson Contouring
Pliers

LD33-150
Slim Arch Forming
Plier, 12 cm
Max: .022"x0.25"
(.51mm x .64mm)
Wires

Clear Brace Aligner Pliers



LD34-200
Hole Punch
Cutter

LD34-201
Tear Drop
Cutter


LD34-202
Horizontal
Dimple

LD34-203
Mini Dot
Dimple

LD34-204
Vertical
Dimple

LD34-205
Dimple Remover

Individual Torquing Pliers Set

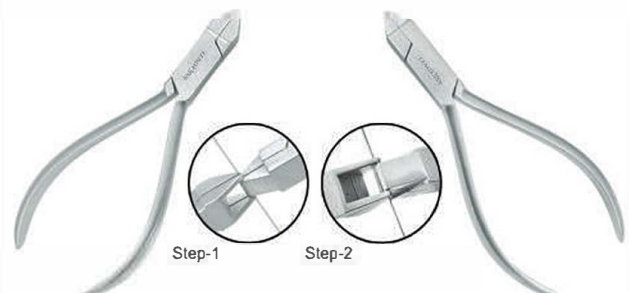

LD33-156
Individual Torquing Pliers set
To make 6mm long torque 13cm
Can bend wire upto. .022" x .025"

Hard Wire Cutters



LD36-101
Hard wire Cutter
Straight



LD36-102
Hard wire Cutter
15° Curved



LD36-103
Hard wire Cutter
Black 15° Curved



LD36-104
Hard wire Cutter

WIRE CUTTER TC GOLD

Tungsten Carbide (Tip Insert)



LDD36-2-16 cm
LDD36-3-18 cm
Wire Cutter
TC Gold
Double Action

Hand Instruments



LD36-105



LD36-106



LD36-107



LD36-108



LD36-109



LD36-110



LD36-111

Uses of Articulator

Primary purpose:

- 1) To hold opposing casts in a predetermined fixed relationship
- 2) To open and close.
- 3) To produce border and inter-border diagnostic sliding motions of the teeth similar to these in the mouth.



LDAR-1001

Technical Datas:

Intercondylar distance 120 mm. / condylar inclination fixed at 30°
Bennet angle fixed at 7° / Incisal table inclination 10°
Distance between mounting plates 85 mm. / Weight 400 gr.
Opening 180°



LDAR-1002
Black



LDAR-1003
Blue



LDAR-1004
Yellow

Impression Trays Dentulous Perforated



Impression Trays Edentulous Perforated



Impression Trays Dentulous Perforated



Impression Trays Edentulous Non-Perforated



Impression Trays Dentulous Perforated



Impression Trays Edentulous Perforated



Impression Trays Dentulous Non-Perforated



Lenox Impression Trays Sizes Chart

Uppers / Superiori Sizes / Misure (mm)						
	U1 XXS	U2 XS	U3 S	U4 M	U5 L	U6 XL
	56	65	68	71	73	78
	44	52	54	58	59	63
Lowers / Inferiori Sizes / Misure (mm)						
	L1 XXS	L2 XS	L3 S	L4 M	L5 L	L6 XL
	56	66	70	73	76	78
	44	50	52	56	58	60

All Lenox Ins Impression Trays are Made by Surgical Grade Stainless Steel



LD37-105
U1P / 12,5cm
For upper left jaw
For lower right jaw

LD37-106
U1P / 12,5cm
For upper right jaw
For lower left jaw



LD37-107
U1P / 12,5cm
Adjustable Perforated



LD37-108
U1P / 12,5cm
Adjustable Solid



LD38-101
Cotton Holder



LD38-102
Spirit Lamp



LD38-103
Endodontic Box



Kidney Trays

With lid available	Without lid available
LD38-104 6"	LD38-108 6"
LD38-105 8"	LD38-109 8"
LD38-106 10"	LD38-110 10"
LD38-107 12"	LD38-111 12"



LD38-112 6" x 8"
LD38-113 8" x 10"
LD38-114 10" x 12"
With Lid

LD38-115 6" x 8"
LD38-116 8" x 10"
LD38-117 10" x 12"
Without Lid



LD38-118
Instruments Box
23x13x5cm




LD38-119
SCALER TRAY
10 x 20 mm



LD38-120
Dental PRF Box



LD39-101
Dental Scaler Clip Cassettes
for 10 Pcs

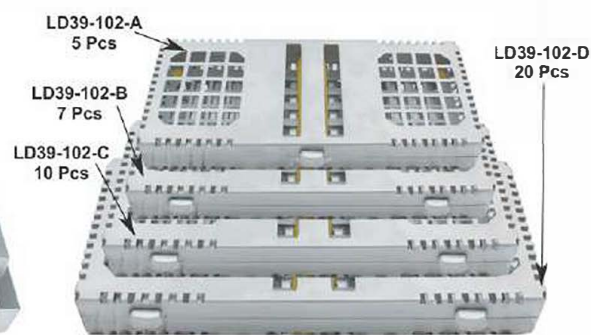

LD39-102-A
5 Pcs
LD39-102-B
7 Pcs
LD39-102-C
10 Pcs

LD39-102-D
20 Pcs

LD39-102


Dental Scalers Cassettes, A, B, C, D



LD39-103
Dental Scaler Clip Cassettes
for 10 Pcs



LD39-104
Dental Scaler Cassettes
for 10 Pcs

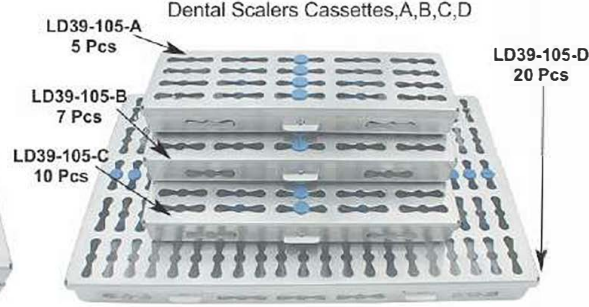

LD39-105-A
5 Pcs
LD39-105-B
7 Pcs
LD39-105-C
10 Pcs

LD39-105-D
20 Pcs

LD39-105


Dental Scalers Cassettes, A, B, C, D



LD39-106
Dental Instruments Cassettes
for 5 Pcs



LD39-107
Dental Scaler Cassettes
for 10 Pcs

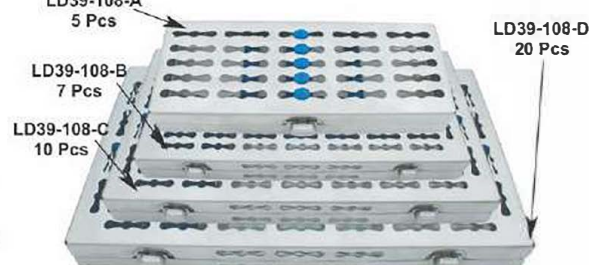

LD39-108-A
5 Pcs
LD39-108-B
7 Pcs
LD39-108-C
10 Pcs

LD39-108-D
20 Pcs

LD39-108

Dental Scalers Cassettes, A, B, C, D

**The Oral Surgery Kit #1001:**

1: 3.0mm Aspiration Tip / 2: Serrated Lucas Surgical Curette, Small
 3: Serrated Lucas Surgical Curette Large / 4: Miller #1 Bone File / 5: Cotton Pliers
 6: Medium Straight Elevator / 7: #77R Back Action Elevator 8: Small Straight Elevator
 9: Lower Universal Apical Forceps / 10: Lower Cowhorn Forceps
 11: Upper Universal Forceps / 12: Lower Universal Forceps / 13: Cone Socket Mirror Handle
 13: Bard-Parker Scalpel Blade Handle / 14: Curved Kelly Haemostat / 15: #5 Cone Socket Mirror
 16: Cottle-Wood Needle Holder / 17: Seldin Flap Retractor / 18: Molt 9 Periosteal Elevator
 19: UNC Peno-Probe / #23 Explorer / 20: Minnesota Retractor / 21: Blumenthal Rongeur
 22: Super-Cut Kelly Scissors / 23: Adson/Brown Tissue Forceps, 2x8 Row Teeth
 24: Extra Large Cassette

**The Soft Tissue Kit #2002****Description**

1: Perio Chisel #13 / 2: Large Cassette / 3: Locking Cotton Pliers
 4: Drop Blade Scalpel Handle / 5: Cone Socket Mirror Handle
 6: #5 Cone Socket Mirror / 7: Molt Curette 2/4 / 8: Cottle-Wood Needle
 Holder / 9: Sugarman 3/4 Perio File / 10: Goldman-Fox #11 Perio Knife
 11: Perio-Probe / #23 Explorer
 12: Younger-Good 7/8 Curette
 13: Goldman Fox Scissors

**Micro Surgery Kit #3003**

1: Perio Chisel #13L / 2: Perio Chisel / 3: Mesial Cassette
 4: Micro Cooley Forceps / 5: Twist Joint Castro-Viejo Needle Holder
 6: Periosteal Elevator / 7: Castro-Viejo Scissors
 8: Molt 9 Periosteal Elevator

**The Implant Site Dev Kit #4004**

1: Small Distal Periosteal TT Coated / 2: Small Mesial Periosteal TT Coated
 3: Serrated Lucas Surgical Curette, Small / 4: 3.5mm Graft Syringe
 5: Large Cassette / 6: Upper Universal Deep Grip Forceps
 7: Lower Universal Deep Grip Forceps / 8: Bone Condenser #12
 9: Cottle-Wood Needle Holder / 10: Mini Buser Periosteal Elevator
 11: Super-Cut Iris Scissors / 12: Suture Forceps, 1.6mm

**LD Hygiene Kit 5005:**

1: Medium Cassette / 2: Cone Socket Mirror Handle
 3: #5 Cone Socket Mirror / 4: UNC Perio-Probe, Color
 5: #204S Scaler / 6: Gracey 11/12 Curette
 7: Gracey 13/14 Curette / 8: Columbia 4R/4L Curette
 9: Gracey 1/2 Curette / 10: H6/7 Scale

**Basic Exam Kit 6006**

1: Small Cassette / 2: Cotton Pliers
 3: Cone Socket Mirror Handle
 4: #5 Cone Socket Mirror
 5: UNC Perio-Probe / #23 Explorer

**Basic Exam Pack: 7007**

1: Regular #5 Mirror & Handle, Solid Handle
 2: Explorer 17/23, Solid Handle
 3: Mariam Cotton Forceps

STEP ONE:

Following clinical use, all instruments must be thoroughly rinsed with water to remove all residue prior to instrument processing.

STEP TWO:

Instruments should be placed into an ultrasonic bath for a period of 20 minutes to ensure complete instrument cleaning.

Note: Ultrasonic Solution should be changed daily or when the solution appears to be visibly soiled.

STEP THREE:

All instruments must be rinsed with water and completely dried before autoclaving. All stainless-steel instruments should be placed in separate sterilization pouches and not to be mixed with any carbide instruments during instrument processing.

NOTE: Only demineralized water should be used in autoclave equipment according to autoclave manufacturer instructions.

- ▶ All stainless-steel instruments must be cleaned and thoroughly dried prior to sterilization.
- ▶ We recommend using only a non-corrosive detergent.
- ▶ Abrasive brushes or materials should not be used to clean instruments.

2-STAINLESS STEEL INSTRUMENTS

We use many different types of stainless steels in the production of instruments dependent upon the design and application of the instrument. All stainless steels have a high nickel and chromium content to maximize corrosion resistance but will corrode and discolor when subjected to high concentrations of certain chemicals.

Do Not expose stainless steel dental instruments to the following chemicals. These chemicals will cause an adverse reaction and may destroy your instruments: Chlorine or Chlorinated products, Household Bleach, Tarter and Stain Remover, Aluminum Chloride, Aqua Regia, Barium Chloride, Bichloride of Mercury, Calcium Chloride, Carbolic Acid, Chlorinated Lime, Citric Acid, Dakin's Solution, Ferric Chloride, Ferrous Chloride, Hydrochloric Acid, Iodine, Lysol®, Mercury Chloride, Mercury Salts, Phenol, Potassium Permanganate, Potassium Thiocyanate, Sodium Hypochlorite (bleach), Stannous Chloride, Sulfuric Acid and Tartaric Acid (Tarter & Stain Remover)

Note: Be sure to only use solutions and chemicals that are compatible with stainless steel hand instruments.

Do Not use sodium chloride solution to clean the instruments. Never exceed temperatures 350o F / 177o C as this will have an adverse effect on the temper of the steel.

Recommended:

Steam sterilization in an autoclave is the benchmark process and is recommended. Sterilization requires a constant temperature of 134°C for 18 to 20 minutes.

3 TITANIUM NITRIDE COATED INSTRUMENTS

Titanium coatings are used to increase the surface hardness of instrument tips and reduce "pullback" when manipulating composite materials. Proper care should be taken to remove any residual composite materials from the blade within 5 minutes after use by wiping the blade with alcohol on a 2x2 gauze. Composite materials or residue can harden on the blade and affect the quality of future restorations. Composite residue left to harden on the blade can not be removed without damaging the surface finish of the composite placement instrument. Titanium Nitride Coated instruments can be cleaned and sterilized using the same recognized acceptable methods as stainless steel instruments

Note: Never use abrasives to clean the surface of any titanium coated composite placement instruments, elevators, periotomes or curettes.

HINGED INSTRUMENTS

All hinge instruments should be allowed to thoroughly dry in an open position after ultrasonic and manual cleaning procedures. We recommend that you only use lubricants such as surgical milk specifically formulated for dental and surgical instruments and follow manufacturer's instructions for applications.

Thank You.

LENOX INS.

LENSUR INSTRUMENTS LTD.

Our Brands



LENOX SURGICAL INSTRUMENTS

Lenox Street #1, Behind Honda Falcon
Daska Road, Sialkot. 51310 Pakistan.
Tel: +92 - 52 - 3241123
Call: + 92 - 320 - 6338182
E-mail: lenoxsurgical@gmail.com
URL: www.lenoxsurgical.com

Authorized Distributors/Representative

